# OMRON 705IT

**Instruction manual**  
**Page 2**

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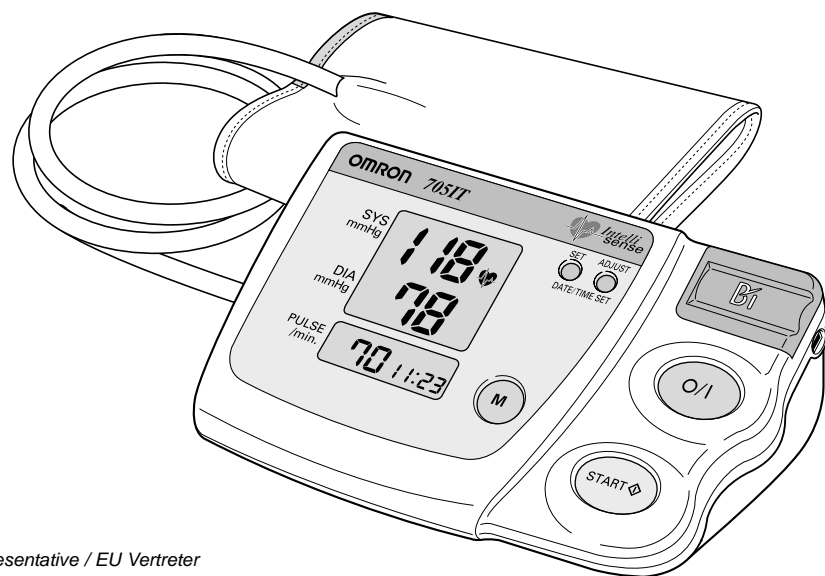
*EU representative / EU Vertreter*  
*Représentant UE / Rappresentante UE*  
*Representante UE / EU Vertegenwoordiger*

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**Fully automatic blood pressure monitor**

705IT Int-0-0

## Contents

Important instructions to help you obtain meaningful readings .....	3
Your blood pressure monitor: the benefits .....	4
Battery installation / replacement .....	5
How to set time and date .....	6
How to apply the cuff .....	7
How to take a reading .....	9
How to use the memory function .....	11
How to use the computer interface .....	12
Hints on taking blood pressure readings .....	16
Failure, causes and rectification .....	17
Maintenance and storage .....	18
Technical data .....	19
OMRON spare parts .....	20
Some useful information about blood pressure .....	21

## Dear Customer

Thank you for purchasing an OMRON product. You can be sure you have made a wise choice. By buying the OMRON 705IT blood pressure monitor, you have chosen a high-quality, innovative instrument for health monitoring. During its development, particular emphasis was placed upon reliability and ease of use. **Before using for the first time, please read through this manual carefully.** If you should still have any questions regarding its use, please contact the OMRON distributor at the given address on the packaging. They will be pleased to help you.

Best wishes for a healthy future.  
Yours sincerely,

OMRON HEALTHCARE

Regular blood pressure monitoring is essential for the prevention, the control and the management of hypertension. It also helps to support the doctors' work.



The World Hypertension League is a world-wide association of experts specialising in hypertension. The World Hypertension League recommends regular blood pressure monitoring by doctors and patients with clinically validated devices.

As the worlds leading manufacturer OMRON provides a range of products for regular blood pressure monitoring.

The World Hypertension League recommends regular blood pressure monitoring by doctors and patients with clinically validated devices





## Important instructions to help you obtain meaningful readings


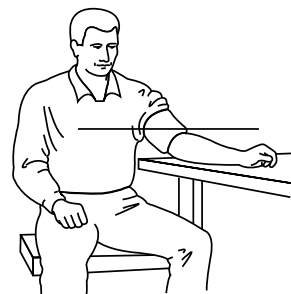
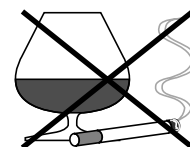
- **Important:** Blood pressure measurement is not suitable in cases of serious arteriosclerosis (hardening of the arteries).
- The pulse display is not suitable for monitoring the frequency of cardiac pacemakers.
- If you suffer from disorders of heart rhythm, known as arrhythmia, you should only use this blood pressure monitor in consultation with your doctor. In certain cases the oscillometric measuring method can produce incorrect readings.
- Pregnant women should only measure their own blood pressure in consultation with their doctor, since the readings may be changed by pregnancy.
- You should avoid eating, drinking (alcohol), smoking and sports before measuring your blood pressure, as this could affect your blood pressure level.
- Don't move, don't speak while measurement is being taken.
- Make yourself comfortable in a chair and relax before each reading.
- We recommend that you monitor your blood pressure twice a day, in the morning after getting up and in the evening after work, or as advised by your doctor.
- **Please remember: Self-measurement is not the same as medical treatment !**  
**You should never change the dose of medicines prescribed by your doctor.**

## Your blood pressure monitor: the benefits

OMRON 705IT is a compact, fully automatic blood pressure monitor, operating on the oscillometric principle. It measures your blood pressure and pulse rate simply and quickly. It contains an intelligent "fuzzy logic" system for controlled inflation, known as "Intellisense". This is an advanced method of oscillometric measurement and does not require pressure pre-setting or re-inflation. OMRON is the world's leading manufacturer of blood pressure monitors with the goal to always fulfill the need for reliable regular monitoring.
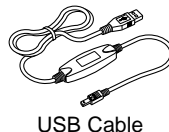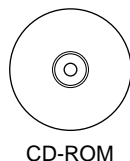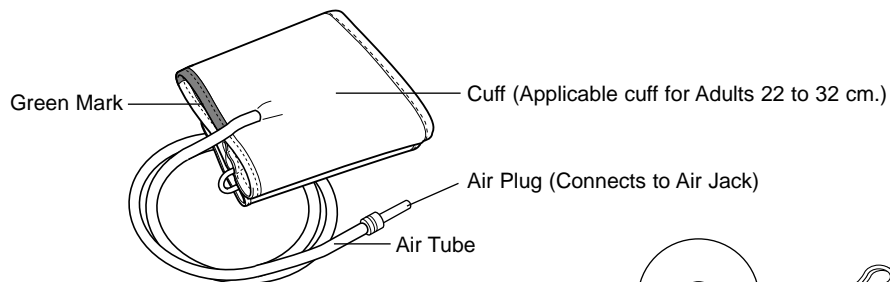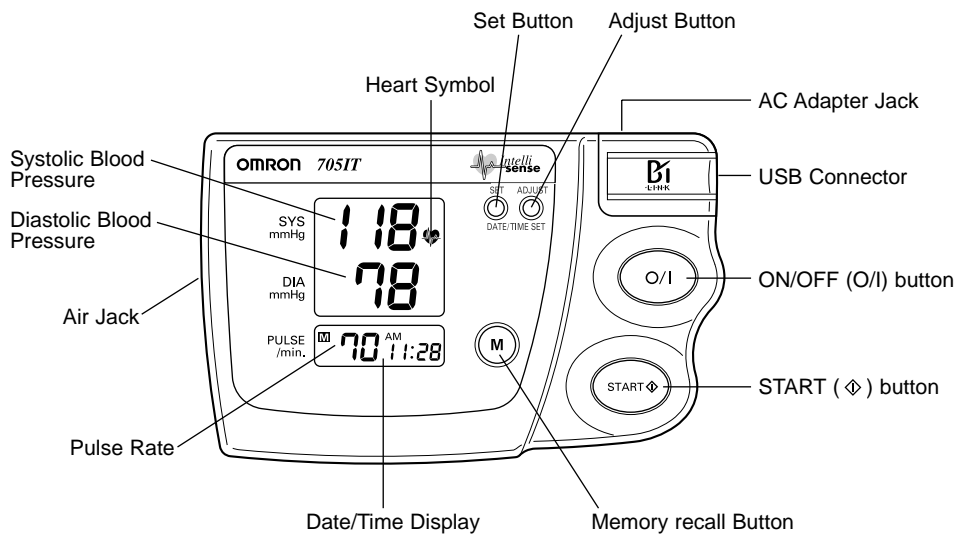
## Battery installation / replacement

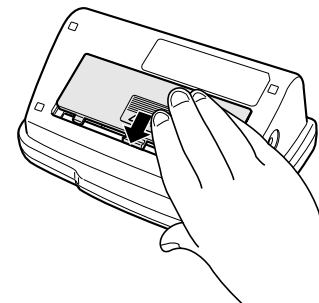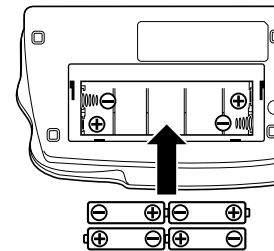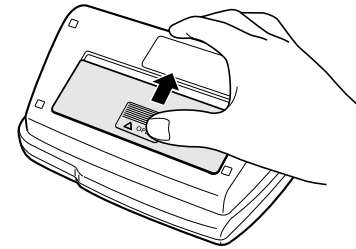- 1 Slide the battery cover off in the direction of the arrow.
- 2 Install or replace 4 "LR6" batteries so the + (positive) and - (negative) polarities match the polarities of the battery compartment as indicated.
- 3 Replace the battery cover.






### NOTES :

If the Low Battery Indicator ( ) appears on the display, replace all four batteries. Long-life alkaline batteries are recommended.

Remove the batteries if the monitor will not be used for an extended period of time.

Replace batteries within 30 seconds. If the batteries remain removed longer than that time, the set date and time will be deleted.


## How to set time and date

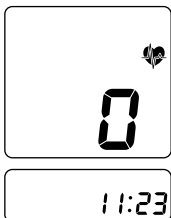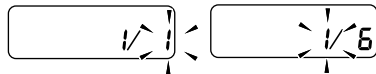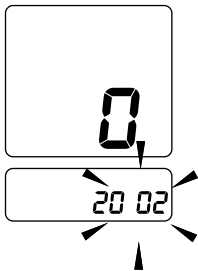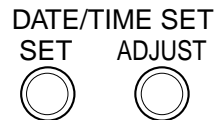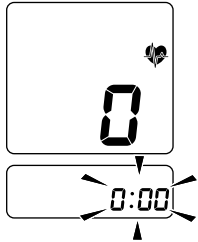- ❶ When the batteries are installed, the display will show 0:00.
- ❷ Press and hold the SET button. The year digits (2002) will flash.
- ❸ Press the ADJUST button to advance the digit(s) one at a time. If you hold down the ADJUST button, the digits will advance rapidly.
- ❹ Press SET button when desired number is on the display to lock setting.
- ❺ Repeat Step 3-4 for month and date.
- ❻ Repeat Step 3-4 for hour and minutes.


## How to apply the cuff

(when wrapping on your left arm)

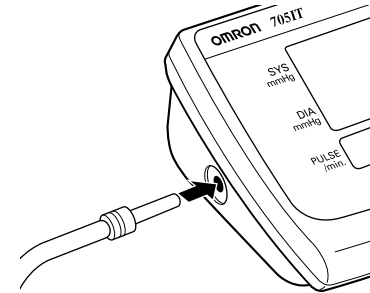
- You can wrap the cuff either on your right or left arm.

- ❶ Insert air tube to air jack (on the left side of the device). The cuff must be airless.
- ❷ Remove tight-fitting clothing from your upper arm.
- ❸ Sit in a chair with your feet flat on the floor and place your arm on a table so that the cuff is at the same level as your heart.
- ❹ Put your arm through the cuff loop. The bottom of the cuff should be approximately 1 to 2 cm above the elbow. The green marker on the cuff should lie over the brachial artery on the inside of the arm. The tube should run down centre of arm approximately even with middle finger.

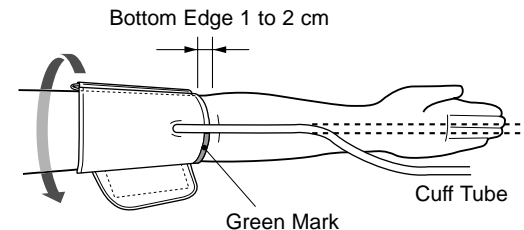
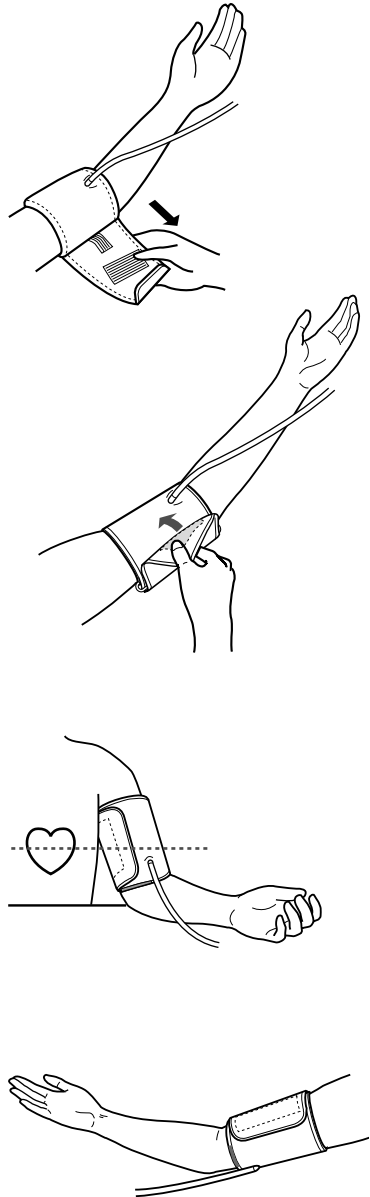
Right hand measurement

- 5 Pull the cuff so that the top and bottom edges are tightened evenly around your arm.
- 6 When the cuff is positioned correctly, close the velcro fastener **FIRMLY**.
- 7 Make certain the cuff fits snugly around your arm. The cuff should make good contact with your skin. **You should be able to fit your index finger between the cuff and your arm easily, so you can pull the cuff off and on.**
- 8 Relax your arm and turn your palm upward.
- 9 Be sure there are no kinks in the air tubing.

## How to take a reading

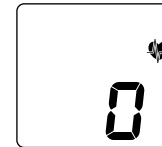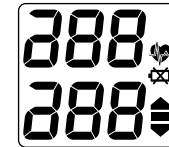- 1 Press the (O/I) button.
  - a) All display symbols appear for approximately one second.
  - b) When the monitor becomes ready to measure, the (♥) symbol appears on the display.
- 2 Press the (START ◊) button and remain still.
 

As the cuff begins to inflate, the monitor automatically determines your ideal inflation level. Because this monitor detects the pulse even during inflation, do not move your arm but remain still until the entire measurement completes.

  - If you want to stop the inflation or measurement, press the (O/I) button. The monitor will stop inflating and start deflating rapidly, then the monitor will turn off.
- 3 Inflation stops automatically and measurement is started.
 

As the cuff slowly deflates, decreasing numbers appear on the display and the (♥) symbol flashes at every heartbeat. In rare circumstances, a higher inflation may be necessary. In those cases, the monitor reinflates the cuff up to 30 mmHg higher than initial inflation and restarts the measurement.

- ④ When the measurement is complete, the cuff completely deflates and your blood pressure and pulse rate are displayed.

**NOTE :**

The monitor automatically stores blood pressure and pulse rate into the memory.

- When 28 sets of readings are stored in memory, the oldest set will be deleted to store a new set.

- ⑤ Press the (O/I) button to turn the monitor off.

**NOTE :**

If you forget to turn the monitor off, it will automatically shut itself off after five minutes.


## Instructions for special conditions:

- If your systolic pressure is known to be more than 220 mmHg, press and hold the (START ◊) button until the monitor inflates 30 to 40 mmHg higher than your suspected systolic pressure.

**NOTE :**

**Do not apply more pressure than necessary.**

**The monitor will not inflate above 300 mmHg.**

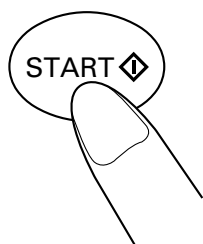
Continue to press button

## How to use the memory function

This monitor has a memory capable of storing 28 sets of readings. Every time you complete the measurement, the monitor automatically stores blood pressure and pulse rate.


- Turn on the monitor by pressing the (O/I) button.
- Press the (M) button to recall stored readings while the (P) symbol is displayed.

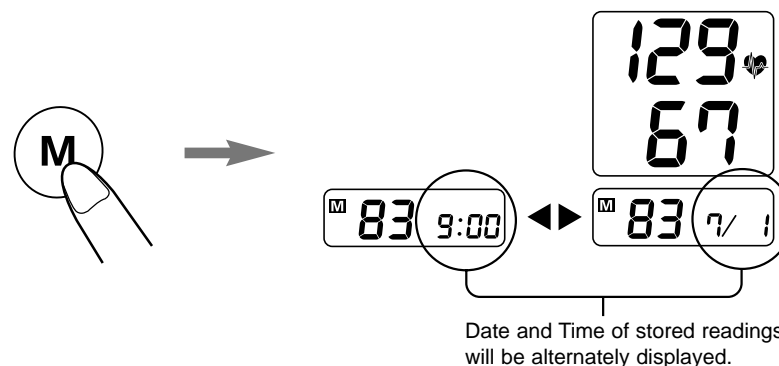- To delete all the stored readings, press the (M) button and the (START ◊) button simultaneously.

**NOTE :**

You can not delete the stored readings partially.
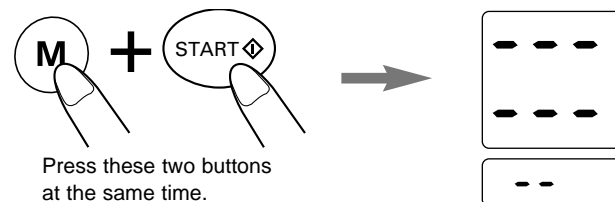
## How to use the computer interface

This blood pressure monitor has a computer interface and is supplied with software to transfer the memory of the blood pressure monitor to your personal computer.

### Required operating environment

The following operating environment is required to operate the communication software as supplied on the CD-ROM.

Personal computer	:	PC/AT compatible with CPU of 80486 or above or with any Pentium processor.
OS	:	Windows 98SE/ME/XP/2000
Memory (RAM)	:	16 MB or more
Hard disk	:	Free area of 10 MB or more is necessary
Display	:	640 x 480 pixels or more
Disk drive	:	CD-ROM drive
Connection port	:	Free USB connection port

### How to install the software

- 1) Turn on your personal computer  
.....Your Windows™ desktop is displayed
- 2) Insert the CD-ROM in the disk drive
- 3) Open the control panel
  - Click the “Start” button.
  - Point the mouse pointer to “Settings” on the menu.
  - Click “Control Panel”
- 4) Open the “Add/Remove Programs” by double clicking the icon
- 5) Execute the installation
  - Click the “install button”
  - Follow the instructions
  - The software is installed on your personal computer

### How to uninstall the software

- 1) Open the control panel
  - Click the “Start” button
  - Choose “Settings”
  - Click “Control Panel”
- 2) Open the “Add/Remove Programs” icon by double clicking
- 3) Delete the application
  - Click “BPFile” from the list
  - Click the Add/Remove button
  - .....The confirmation screen is displayed
  - Click “Yes”
  - .....“BPFile” is deleted
  - Click “OK”

### How to connect your 705IT blood pressure monitor to your PC

Connect your 705IT blood pressure monitor with your PC using the included USB cable. The mini USB connector (small connector) should be connected to your 705IT blood pressure monitor and the normal USB connector (large connector) should be plugged into your personal computer.

Please ensure you connect the USB cable the correct way to your 705IT blood pressure monitor, with the [O] symbol facing upwards.

## Hints on taking blood pressure readings

- Before measuring, stretch your back and sit up straight, breath normally 5-6 times, relax your shoulders, arms and entire body.
- Slightly bend and support your elbow (e.g. on a table).
- Do not move, talk or touch the device during measurement.
- Stress influences your blood pressure, so try to relax before taking a measurement.
- To compare results, measurements must always be taken from the same arm.
- Do not place the cuff over thick clothes and do not roll up your sleeve if it is too tight.
- Align the cuff with the height of your heart (breast height). If necessary use a pillow or cushion.
- Coldness raises your blood pressure. It is recommended to measure at room temperature (+20°C).
- Relax before measurement. If necessary go to the bathroom before measurement.
- Do not measure directly after bathing.
- If you wish to interrupt measurement prematurely, press the (O/I) button. Inflation is stopped and the air is automatically released from the cuff.
- Blood pressure measurements taken in a moving vehicle may be influenced by vibration.
- If an error has occurred during measurement **E** is displayed. Please refer to the next chapter.
- The monitor should be switched off after measurement to avoid running down the batteries unnecessarily. If you forget to do so, the monitor switches itself off automatically after 5 minutes.

## Failure, causes and rectification

Failure and their possible causes	Rectification
<p>Incorrect measurement is indicated by <b>E</b></p> <ol style="list-style-type: none"> <li>1. A correct reading could not be obtained because measurement was disturbed by movement of the body.</li> <li>2. The cuff is not fitted correctly.</li> <li>3. Your clothing has constricted your blood flow.</li> <li>4. There is still air in the cuff when the monitor is switched off.</li> </ol>	<ol style="list-style-type: none"> <li>1. Repeat the measurement keeping perfectly still. Do not move your arm and do not speak.</li> <li>2. Check that the cuff is correctly fitted. Then repeat the measurement.</li> <li>3. Remove the item of clothing which caused the constriction.</li> <li>4. The unit may be defective. Please have it examined and, if necessary, repaired by OMRON.</li> </ol>
<p><b>The (START  ) button was pressed before the (  ) symbol was displayed.</b></p>	<p>Wait for the (  ) symbol indicating readiness before you press the (START  ) button.</p>
<p><b>The blood pressure values displayed are extremely low or high, or they are implausible.</b></p>	<p>Refer to the instructions and then repeat the measurement.</p>
<p><b>The display does not light up when the (START  ) button is pressed.</b></p> <ol style="list-style-type: none"> <li>1. The batteries are exhausted.</li> <li>2. The +/- poles of the batteries have been reversed.</li> <li>3. The battery contacts are dirty.</li> </ol> <p> The batteries are weak or exhausted.</p>	<ol style="list-style-type: none"> <li>1. Check the batteries and, if necessary, fit 4 new batteries.</li> <li>2. Reinsert the batteries with the +/- poles the right way round.</li> <li>3. Clean the battery contacts with a dry cloth.</li> </ol> <p>Fit new batteries.</p>
<p><b>The cuff pressure does not rise although the pump motor can be heard.</b></p> <p><b></b> problem with mem. function</p>	<p>Check that the air tube is properly connected to the monitor. Push the tube connector firmly into the socket.</p> <p>Contact your Omron retailer or distributor.</p>

## Maintenance and storage

- For cleaning the monitor only use a soft, lightly moistened cloth.  
**Do not use petroleum spirits, thinners or similar solvents !**
- Stains on the cuff can be carefully removed using a moist cloth and soapy water.  
**Do not wash the cuff or make it wet !**
- Keep your blood pressure monitor in the storage case to protect it from dust and moisture.
- Protect your monitor against vibrations and do not leave it in a place where temperatures are extremely low (below -20°C) or extremely high (above 60°C) or with high humidity (+ 85%).
- Do not use your monitor at very low temperatures (below 10°C) or very high temperatures (above 40°C).
- Do not fold the cuff or air tube together too tightly and do not crease them.
- Do not carry out repairs of any kind yourself. If a defect occurs, consult the OMRON distributor or Customer Services or seek the advice of your surgical appliance stockist or pharmacist.
- The accuracy of this blood pressure monitor is designed to last a long time.
- In order to increase the lifetime of the batteries take them out when you store the device for a longer period than 2 months.
- Electro-magnetic fields can influence the measurement. Do not use a portable phone near this unit.


## Technical data

Model:	OMRON 705IT
Display:	LCD digital display
Blood Pressure	Pressure: 0 mm Hg to 299 mm Hg
Measurement Range:	Pulse: 40 – 180/min.
Accuracy:	Pressure: ±3 mm Hg Pulse: ±5% of reading
Memory:	28 measurements
Inflation:	Fuzzy-Logic controlled by electric pump
Deflation:	Automatic pressure release valve
Pressure Detection:	Capacitive pressure sensor
Measurement Method:	Oscillometric method
Power Source:	4 Alkaline batteries 1.5V (Type LR6) or AC/DC adapter (optional, 6V = 4W)
Battery Life:	Capacity of new batteries is approx. 300 measurements
Operating Temperature:	10°C to 40°C
Storage Temperature:	-20°C to 60°C
Humidity:	15 to 90% RH maximum
Storage humidity:	10 to 95% RH maximum
Console Weight:	Approximately 380g without batteries
Outer Dimensions:	Approximately 115 (l) mm x 177 (w) mm x 71 (h) mm
Cuff Dimensions:	Approximately 140 mm x 480 mm (arm circumference 22 – 32 cm)
Accessories:	Medium cuff, CD-ROM, USB cable, instruction manual, storage case, battery set

This device fulfils the provisions of the EC directive 93/42/EEC (Medical Device Directive). This blood pressure monitor is designed according to the European standard EN1060, Non-invasive sphygmomanometers Part 1: General requirements and Part 3: Supplementary requirements for electromechanical blood pressure measuring systems.

**NOTE:** Subject to technical modification without prior notice

- Disposal of this product and used batteries should be carried out in accordance with the national regulations for the disposal of electronic products.

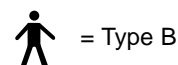
CE 0197

## OMRON spare parts

Printer

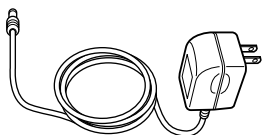AC adapter

Small cuff  
Arm circumference 17 – 22 cm

Medium cuff  
Arm circumference 22 – 32 cm

Extra large cuff  
Arm circumference 32 – 42 cm

## Calibration

The accuracy of this blood pressure monitor has been carefully tested and is designed for a long service life. It is generally recommended to have the monitor inspected every two years to ensure correct functioning and accuracy. Please consult your authorised dealer or the OMRON Customer Service at the address given on the packaging/attached literature.

OMRON devices are especially designed for regular blood pressure monitoring. OMRON therefore has direct contact with medical specialists and takes advice, on the design and features of its blood pressure monitors.

## Spare parts & accessories

Due to high company quality standards, OMRON considers the main unit as a non-serviceable part because of the necessity of proper calibration after replacement of high-tech components.

## Some useful information about blood pressure

### What is Blood Pressure ?

Blood pressure is a measure of the force of blood flowing against the walls of the arteries. Arterial blood pressure is constantly changing during the course of the heart's cycle. The highest pressure in the cycle is called the **SYSTOLIC BLOOD PRESSURE**; the lowest is the **DIASTOLIC BLOOD PRESSURE**. Both pressure readings, the **SYSTOLIC** and **DIASTOLIC**, are necessary to enable a doctor to evaluate the status of a patient's blood pressure.

Many factors such as physical activity, anxiety, or the time of day, can influence your blood pressure. Blood pressure is typically low in the morning and increases from afternoon to evening. It is lower in the summer and higher in the winter.

Blood pressure is measured in millimetres of mercury (mmHg) and measurements are written with the systolic pressure before the diastolic e.g. A blood pressure written as 140/90, is referred to as 140 over 90 mmHg.

### High Blood Pressure

The discovery that you have high blood pressure is more often than not a chance finding. People with high blood pressure usually feel well unless the blood pressure has been high for some time, and complications have occurred. Such complications may effect the heart, brain and other important organs. They may bring on ill health and affect your enjoyment of life. It is because of the complications, rather than the high blood pressure itself, that this condition is taken so seriously today.


High blood pressure is sometimes referred to as raised blood pressure, elevated blood pressure, or hypertension. Most people with this condition require regular long-term drug treatment, usually in the form of tablets. Some can lower their blood pressure sufficiently by other means, such as changing diet and life-style.

People who have high blood pressure, and who also understand their condition, have an advantage. They are more likely to attend for regular check-ups, to cooperate with advised diets and changes in lifestyle, and to take medicines correctly. This immediately leads to better health. It is hoped that, by a better understanding of high blood pressure and how it is treated, you will take an active and informed part in your health care, and so enjoy a better quality of life.

## What Causes High Blood Pressure ?

Blood pressure is only classed as high if it doesn't go down when you rest. Permanently raised blood pressure can be caused by several factors such as hardening of the arteries, smoking or drinking too much alcohol. It can also be caused by cholesterol, a type of fat, building up inside your blood vessels and making them narrower. Being under stress can add to the problem.

High blood pressure is also known to run in families hence if you find you have a high blood pressure it is a good idea to test other closely related members of your family.
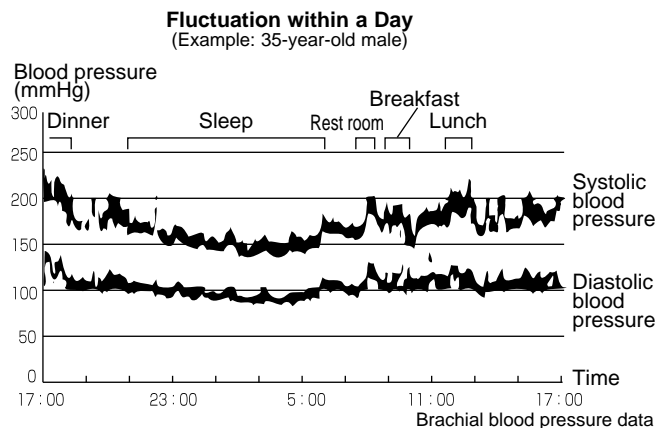
## Why is it a Good Thing To Measure Blood Pressure at Home ?

Having your blood pressure measured by a doctor can cause anxiety which is a cause itself of high blood pressure. As a variety of conditions affect blood pressure, a single measurement may not be sufficient for an accurate diagnosis.


The blood pressure measured first thing in the morning after getting up, before eating and while at rest, is known as the fundamental blood pressure. In practice it is rather difficult to record the fundamental blood pressure, but it is important to get readings in an environment that is as close as possible to this.

## Classification of Blood Pressure by the World Health Organization

The World Health Organization (WHO) and the International Society of Hypertension (ISH) developed the Blood Pressure Classification shown in this figure.

This classification is based on the blood pressure values measured on people in a sitting position in outpatient departments of hospitals.

*\*There is no universally accepted definition of hypotension. However, those having the systolic pressure below 100 mmHg are assumed as hypotensive.*

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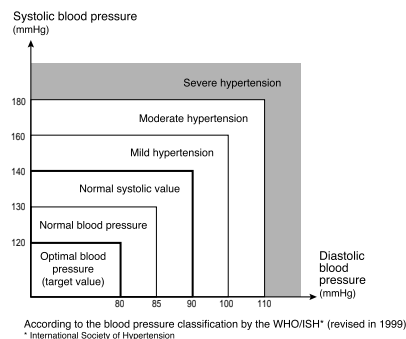## Table of Contents

---

Before measuring your Blood Pressure	3
Notes on safety	4
Names and functions of the parts	5
How to insert the batteries	6
How to set the Date and Time	7
How to measure Blood Pressure correctly	7
How to obtain a correct posture	9
How to measure Blood Pressure	10
How to use the Memory function	11
How to use the Graph function	12
How to change the settings	13
Auto memory	14
Position sensor	15
Left or right wrist measurement	15
Alarm	15
Time	16
Delete data	16
Reset settings	16
What is Blood Pressure?	17
Blood pressure at the upper arm and wrist	18
Maintenance and storage	19
Error indication	20
Troubleshooting	21
Spare parts and accessories	22
Specifications	22

---

Dear Customer,

Thank you for purchasing this OMRON product. OMRON products are known for their reliability and advanced technology. By buying this OMRON R7 wrist blood pressure monitor with position sensor you have chosen for a high quality and innovative instrument. We hope that you will be fully satisfied with your new blood pressure monitor. Please take your time to read this instruction manual carefully before using it for the first time. If you still have questions please contact your OMRON dealer or the distributor as mentioned on the packaging. They will be pleased to help you.

We wish you a healthy future!

## Before measuring your Blood Pressure

### Notes on Safety Please read thoroughly before using this monitor.

- The warning signs and the sample icons shown here are listed for you to use the product safely and correctly as well as to prevent the risk of damage to you and others.
- The icons and meanings are shown as follows.

Warning sign	Contents
 <b>Caution</b>	Indicates matters in which bodily harm or material damage* may arise as a result of incorrect handling.

\* Material damage refers to a wide range of damage involving your house, household goods, domestic animals, and pets.

Examples of signs	
	The  icon indicates prohibitions (what you cannot do). Matters involving actual prohibitions are indicated by text or pictures in or near  . The left icon refers to "prohibition to disassemble".
	The  icon indicates something compulsory (always follow). Matters involving actual compulsory actions are indicated by text or pictures in or near  . The left icon refers to "general compulsion".

 <b>Caution</b>	
Self-diagnosis of measured results and treatment are dangerous. Please follow the instructions of your doctor. <ul style="list-style-type: none"> <li>• Self-diagnosis may worsen the disease.</li> <li>• When the arterial sclerosis proceeds due to diabetics, hyperlipemia, or hypertension, it may cause peripheral circulatory disturbance. For these cases, blood pressure may differ greatly between measurement at the wrist and the brachial.</li> </ul>	
Do not use this monitor for the infant or the person who cannot express one's intention. <ul style="list-style-type: none"> <li>• They may suffer accident or trouble.</li> <li>Do not use this monitor for purposes other than measuring blood pressure.</li> <li>• You may suffer accident or trouble.</li> <li>Do not use portable phones near this monitor.</li> <li>• The monitor may operate erroneously.</li> </ul>	
Do not disassemble or modify the main unit or the wrist cuff of the blood pressure monitor. <ul style="list-style-type: none"> <li>• Correct measurement will be impossible.</li> </ul>	

## Notes on safety

---

### Important instructions for obtaining meaningful readings

- Avoid eating and drinking (alcohol), smoking, sport and taking medicines before measuring your blood pressure as this could affect your blood pressure level.
- Before you start sit comfortably on a chair and relax. It is recommended to check your blood pressure at least twice a day, in the morning after getting up and in the evening after work, or as advised by your doctor.
- Do not move and do not speak while the measurement is being taken. Avoid measuring your blood pressure while you are in a moving vehicle.
- Allow an interval of at least three minutes between readings. Otherwise congestion of blood might occur and distort the readings. Always remove the cuff between readings.
- Bare your wrist and take off your watch and jewellery before applying the cuff. Do not constrict the flow of your blood if you push up your sleeve.

### Cautions

Blood pressure measurement is not suitable in case of serious arteriosclerosis (hardening of the arteries).

- The pulse rate indicator is not suitable for monitoring the frequency of cardiac pacemakers.
- If you suffer from disorders of heart rhythm, known as arrhythmia, you should only use this blood pressure monitor in consultation with your doctor. In certain cases oscillometric measurement methods can produce incorrect readings.
- Pregnant women should only measure their own blood pressure in consultation with their doctor since the readings may be changed by the pregnancy.

***Please remember: Self-measurement is not the same as medical treatment! You should never change the dose of medicines as prescribed by your doctor!***

### Before using your OMRON R7 for the first time

Your blood pressure monitor: the benefits

The OMRON R7 is a compact fully automatic wrist blood pressure monitor. Operating on the oscillometric principle it measures your blood pressure and pulse rate comfortably and quickly from the wrist. It contains an intelligent system for “fuzzy” controlled inflation, known as “Intellisense”. This is an advanced method of oscillometric measurement. No pre-setting of re-inflation is necessary.

The OMRON R7 has a Position Sensor. The Position Sensor is able to determine the optimal height of your wrist (heart level) for the blood pressure measurement. Measuring at the heart level is extremely important since measuring at a higher or lower position will give wrong results.

The OMRON R7 has a large DOT matrix display enabling you to enlarge the digits and showing you graphs. Together with the built in alarm function for same time measurements or pill intake your OMRON R7 is the optimal device to manage your blood pressure.

Main unit


Batteries for monitoring

(Two LR03 alkaline batteries)

Storage case

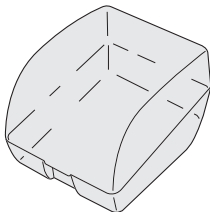Instruction Manual

Blood pressure pass

Guarantee card

## Names and functions of the parts

Main unit

**Systolic Blood Pressure Display**

**Diastolic Blood Pressure Display**

**Date/Time Display**


**Mark**

Displayed when a stored measurement is shown.

**Pulse Display**

**◀▶ Button**

Push ▶ to move the displayed item forward. Push ◀ to move the displayed item backward.

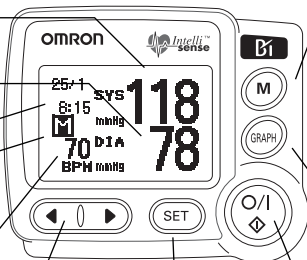**MEMORY RECALL button**

Push this button to see the measured values stored in the monitor.

**GRAPH CALL button**

Push this button to see the graph.

**ON/OFF button**

Push this button to measure your blood pressure or to switch the monitor off.

**SETTING button**

Push this button to set various items on the monitor such as time and date.

## How to insert the batteries




Use this monitor after inserting batteries in the main unit.

1. Open the battery compartment.

2. Insert the batteries (Type LR03 / AAA) as shown in the figure.

Ensure the polarity (+/-) of the batteries is correct.


3. Close the battery cover.

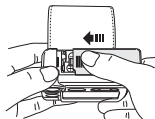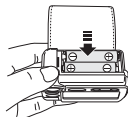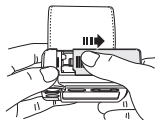With new batteries (long life alkaline batteries) you are able to measure app. 300 times under normal conditions. Please switch off the monitor after using it. If you forget to switch the monitor off it will switch off itself after app. 2 minutes. Please take the batteries out of the monitor if you intend not to use the monitor for a long period (app. 2 months).

## How to set the Date and Time

Before using your OMRON R7 for the first time, you have to set the date and time. After turning the monitor on by pushing the on/off button the initial screen (see picture) will be available for a few seconds. After the initial screen you will proceed automatically to the TIME SET screen. (For instruction how to set the date and time refer also to page 16 in Chapter How to change the settings.)

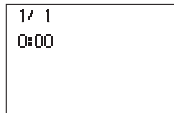Initial Screen

After reaching the time screen you will see the Gregorian calendar (year/day/month). The year will be flashing.

1. By using the ( ◀ ▶ ) buttons you can set the current year.
2. By pushing the SET button you will lock the year and proceed to the next item.

Time Screen

For setting the month, date and time repeat step 1 and 2. After setting the minutes you will proceed to the DELETE DATA screen. Please turn off your monitor by pushing the on/off button. Your OMRON R7 is ready for use now.

## How to measure Blood Pressure correctly


### How to wrap the cuff

In order to measure your blood pressure correctly it is needed to wrap the cuff around your wrist in a correct way. If not, correct measurements are not possible. Please follow the coming instructions (left wrist):



- 1) Make sure that your wrist is bare and will be wrapped over your skin.
- 2) Make sure that your clothing is not constraining the flow of your blood.
- 3) Make sure that the wrist (left or right) to measure is selected from the settings menu (see also page 15). The left wrist is preset.

- 4) Place the cuff over your left wrist with your thumb facing upwards. Leave app 10-15 mm between the cuff and the bottom of the palm of your hand as shown in the figures.

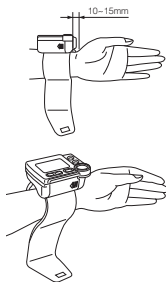- 5) Pull the wrist cuff upward and close the cuff firmly with the Velcro tape as shown.

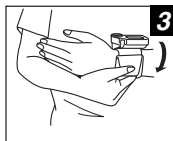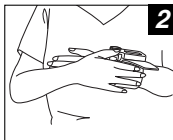You can use your OMRON R7 also on your right hand after selecting right wrist measurement in the SETTING MODE.

People with a slender wrist (wrist circumference of less than 145mm) may feel loose because the cuff may not be wrapped firmly. Please do not force the cuff it will not affect the measurement.

## How to obtain a correct posture




Achieving a correct posture before measuring your blood pressure is important since an incorrect posture might give you incorrect readings. To achieve a correct posture please follow the instructions (see pictures):

- 1) Sit down and stretch your back.
- 2) Relax yourself by relieving the stress from your shoulders and arms.
- 3) Breath several times deeply and return to natural respiration.
- 4) Do not sustain your wrist with your other hand and avoid stress in your hand by making firm grips or large extensions.

### Correct Posture

### Incorrect Posture
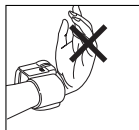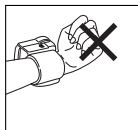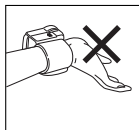
### Measuring your blood pressure and pulse

If you have achieved the correct posture you can start the measurement. When you measure for the first time the Position Sensor and the Signal will be ON and measuring at the left wrist will be pre set.

Push the on/off button to activate your OMRON R7 and the Position Sensor Screen will be displayed.


## How to measure Blood Pressure

Sustain your elbow with your other hand or a table to avoid movement.


• Adjust the height of your wrist by using your elbow as a fulcrum. When the ( ◀ ▶ ) will reach the ( ♥ ) sign you will hear a beeping sound indicating that your blood pressure monitor is at the correct height (heart height). The ( ♥ ) sign will change in the ( 📶 ) sign and the measurement will start.

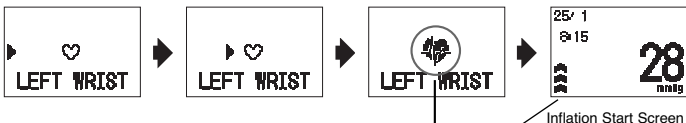**Do not move until the measurement results are displayed!**


The screens (see picture) will be shown during the measurement.

When the measurement is finishes, the deflation mark ( ▼ ) flashed and the cuff deflates.
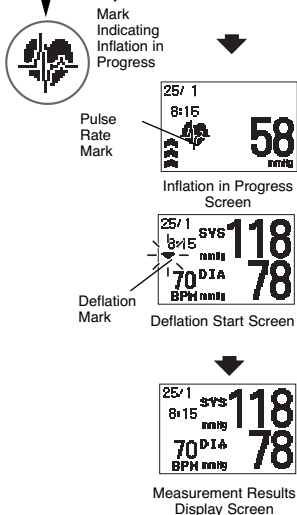
## How to enlarge the reading

The results will be shown on your display. By pushing the ( ◀ ▶ ) buttons you can enlarge the readings.

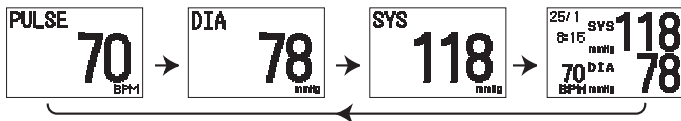
## How to use the Memory function


Your OMRON R7 blood pressure monitor is able to store maximum 90 measurements including date and time. When your blood pressure monitor is holding 90 measurements the oldest will be deleted to store a new measurement.

Press the (M) memory button when your blood pressure monitor is switched off. The latest measurement will be displayed. Each time you push ( ▶ ) button an earlier measurement will be displayed. By pushing the ( ◀ ) button a later reading is displayed.

By pushing the GRAPH button the selected measurement is displayed in a weekly graph (see also How to use the graph function). Push the memory button to return to the memory display again.

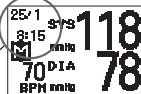Date and Time of Measurement
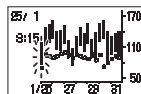
## How to use the Graph function

This monitor has the possibility to display the stored measurements as a graph. In order to get a clear history of your blood pressure the measurements are put in time zones:

- Morning: Measurements made between 0:00 and 9:59
- Daytime: Measurements made between 10:00 and 17:59
- Evening: Measurements made between 18:00 and 23:59

The last measurement in each time zone will be displayed. For example if you measure 3 times in the same time zone on the same day only the third reading will be displayed. On the other hand all readings will be stored in the memory (See also page 11).


There are three different graphs available:

- All data graph: displays all measurements per week
  - Morning data graph: displays all morning measurements per week
  - Evening data graph: displays all evening measurements per week
- No graph display is provided for daytime data only.

To display the graph push the GRAPH CALL button.


All Data Graph is displayed.

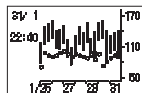All Data Graph

Push the GRAPH CALL button again to display the Morning Data Graph.

Morning  
Data Mark

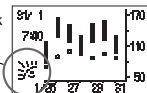Morning Data graph

Push the GRAPH CALL button again to display the Evening Data Graph.

Night  
Data Mark


Night Data Graph


Each time you push the GRAPH CALL button, All, Morning, and Evening data are displayed in that order repeatedly.

## Explanation of the graph display

When the GRAPH button is pushed the latest measurement will appear in a flashing graph display. In the example the measurement is made on January 31 at 22:40. The graph shows 140 mmHg for systolic and 95 mmHg for diastolic (top and bottom of vertical bar) and a pulse rate of 90 beats a minute.

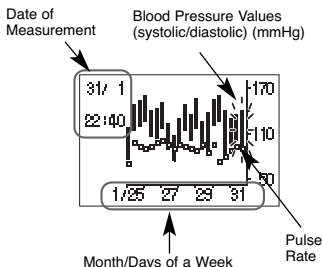When you push the memory button in the graph display the details of the measurement will appear. Pushing the graph button again will bring you back the graph display. You can use ( ◀ ▶ ) the arrows in the graph and memory display to leaf through the measurements.

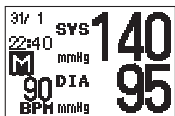
Detail Screen of Data in Flashing Display


## The systolic blood pressure is higher than 170mmHg

When a systolic blood pressure above 170 mmHg is measured twice the ( ▲ ) symbol is displayed in the graph.

When a systolic blood pressure above 170 mmHg is measured three times the ( ▲ ) symbol disappears and the screen automatically changes to the “hypertension display”. The maximum displayed blood pressure value (y-axis) will change from 170 mmHg to 210 mmHg.

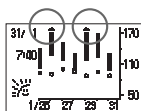Display screen where the systolic blood pressure values higher than 170 mmHg are measured twice

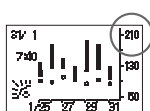Graph Display Screen of Hypertension Type

## How to change the settings


The settings of your OMRON R7 blood pressure monitor can be changed according to your own preferences and comfort.

To enter the setting menu please press the SET button when your blood pressure monitor is switched off. By using ( ◀ ▶ ) you can select your preference. By pressing the SET button again you lock your preference and proceed to the next item as shown in the figure. The on/off button can be used to switch the monitor off again.


• If you push the **SETTING** button while the power is turned off, the Setting Screen is displayed in the order shown in the picture.

• You can select the setting or set the date and time with the buttons.


• To finish the setting, push the on/off button.

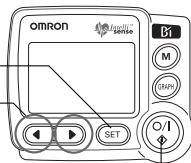
Automatic Memory Setting  
(Refer to Page 14.)

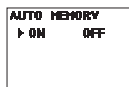
Position Sensor Setting  
(Refer to Page 15.)

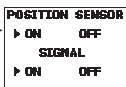
Measuring (Wrist) Setting  
(Refer to Page 15.)

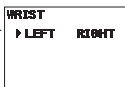
Alarm Setting  
(Refer to Page 15.)

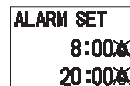
Reset Settings  
(Refer to Page 16.)

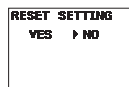All Memory Erasure  
(Refer to Page 16.)

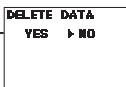Time Setting  
(Refer to Page 16.)
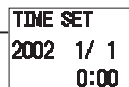
## Auto memory

You can either enable (ON) or disable (OFF) the memory function of your blood pressure. When the auto memory function is switched off the measurement will not be stored. When using your OMRON R7 monitor with more than one person it is recommended to disable the auto memory function. The auto memory function mode is ON when you purchased your OMRON R7 (factory preset).


Automatic Memory  
Screen

Automatic Memory  
OFF Screen

## Position sensor

The position sensor of your OMRON R7 can be switched ON or OFF. In the factory preset the position sensor is enabled. If the position sensor is switched ON you can choose to have the optimal wrist height supported by an acoustic sound (SIGNAL ON or OFF).

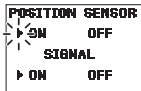Position Sensor Screen


Position Sensor OFF Screen

Position Sensor ON and SIGNAL OFF Screen

## Left or right wrist measurement

Normally blood pressure is measured on the left wrist (factory preset). With your OMRON R7 you can choose to use the monitor either on your left or right wrist. It is recommended to use the same wrist when for optimal comparisons of your blood pressure fluctuations.

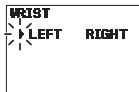Measuring (Wrist) Screen

Measuring (Wrist) Right Wrist Screen

## Alarm


Your blood pressure monitor has the possibility to set an alarm twice a day. The alarm function can be used to monitor your blood pressure at the same time every day or as a reminder for pill intake.

The alarm will sound for 2 minutes. By pushing the on/off button the alarm will stop. The alarm will not sound when the blood pressure monitor is used (e.g. for measuring or graph display) at the moment the alarm is set.

You can set the requested hour by pushing ( ◀ ▶ ) when flashing. Press SET to lock the hour and proceed with the minutes. By pushing the ( ◀ ▶ ) you can set the minutes. Press SET to confirm and to proceed. Push either ( ◀ ▶ ) to enable ( 🔔 ) or disable ( 🔕 ) the alarm. Press SET to lock the settings.

Alarm Screen

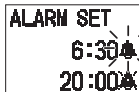Screen of indicating that the alarm is sounding

## Time


You can change the time and date setting of your blood pressure monitor. When reaching the time screen you will see the Gregorian calendar. The year will be flashing.

1. By using the ( ◀ ▶ ) buttons you can set the current year.
2. By pushing the SET button you will lock the year and proceed to the next setting.

For setting the month, date and time repeat step 1 and 2. After setting the minutes you will proceed to the DELETE DATA screen.

Time Screen

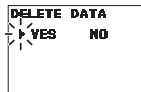
All Memory Erasure  
YES Screen

## Delete data


In this screen you are able to delete all stored measurement from you OMRON R7. Please be careful that deleted data cannot be restored. The factory preset is NO. If you choose YES the monitor will ask for a confirmation "ARE YOU SURE?". To cancel select NO to confirm DELETE DATA select YES.

After deleting all stored measurements the "Data Deleted Screen" will be shown. All the stored measurements are deleted and the monitor will be turned off.

If the stored measurements are not deleted you will proceed to the RESET SETTINGS screen.

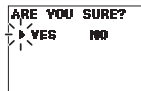Are You Sure YES  
Screen

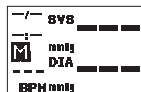
Data Deleted (Blank)  
Screen


## Reset settings

In this menu you can choose to reset all the settings to the "factory preset" settings. The factory preset for this menu is NO. If you choose YES by using the ( ◀ ) button the monitor will ask you for a confirmation "ARE YOU SURE?" After confirmation by selecting YES all the settings will be changed to the factory mode.

The monitor will return to the AUTO MEMORY SCREEN. You can turn off your monitor by pushing the on/off button.

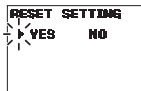Reset Settings YES  
Screen


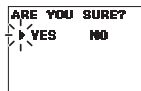
Are You Sure YES  
Screen

## What is Blood Pressure?

Blood is pumped to artery by the beat (contraction and expansion) of a pump called the heart. The pressure to pump out blood to artery is called the "systolic blood pressure." The pressure when the heart expands as the blood having circulated the whole body returns is called the "diastolic blood pressure."

### Health and blood pressure


Blood pressure fluctuates being affected by everyday life and age. Hypertension is a sign that can cause dangerous diseases such as cerebral apoplexy or myocardial infarction to occur easily. By understanding the state of your blood pressure, review your health and lifestyle.
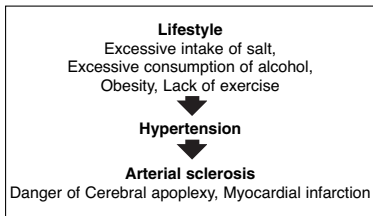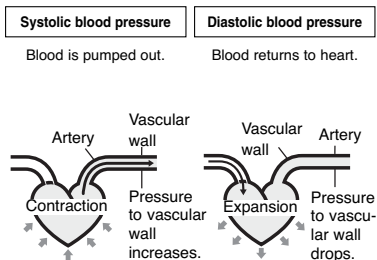
### Hypertension guide

The World Health Organization (WHO) and the International Society of Hypertension (ISH) established the Blood Pressure Classification. Use this Classification as the guide to know the state of your blood pressure. However, if you are concerned with your blood pressure, be sure to follow the instructions of your doctor.

There is no universally accepted definition of hypotension. Generally speaking, those with the systolic pressure below 100 mmHg are considered to be hypotension.

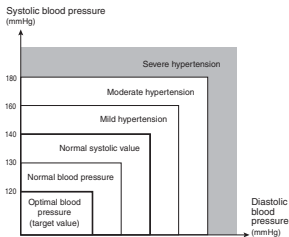According to the blood pressure classification by the WHO/ISH\* (revised in 1999)  
\* International Society of Hypertension

This classification is made based on the blood pressure values measured at upper arm by sitting on the chair in the hospital regardless of the age.


## Blood pressure at the upper arm and wrist

### Wrist Blood Pressure Monitor

Blood pressure is usually measured at brachial artery.

As the wrist arterial vessel is connected to the brachial vessel, wrist blood pressure values are close to brachial ones and reflect the change of brachial blood pressures. However in case there is a circulatory disorder such as arterial infarction, fairly large difference is observed between the brachial and the wrist blood pressure values.

After consulting with your doctor, use the wrist blood pressure together with the brachial pressure for your health management.
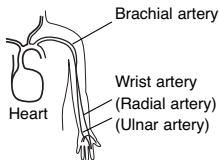
### Diseases influencing blood pressure values

Diabetes, hyperlipemia, or hypertension accelerates arterial sclerosis, which causes dangerous disease such as cerebral apoplexy or myocardial infarction or induces arterial stenosis or peripheral circulatory disorder.

In these cases, the wrist and brachial blood pressure values may differ greatly. (Even in healthy people, a difference of 20 mmHg may occur if the measuring conditions are not proper. Therefore do not try to diagnose the measurement values by yourself, but be sure to follow the instructions of your doctor.)

In addition, difference of wrist and brachial blood pressure values may somewhat get affected by the physiological conditions at the time of measurement. However, the blood pressure values measured at wrist and upper arm usually fluctuate in the same manner. Therefore, you can check the fluctuation tendency of blood pressure by measuring the wrist blood pressure.

### Differences between blood pressure measured at home and at the doctor's office

Normally blood pressure measured at home tends to be lower than when measured at the doctor's office. This can be as much as 20 mmHg due to that you feel more relaxed at home and feel slightly tense when visiting the doctor. To manage your blood pressure it is important to monitor and register your blood pressure at home.

Since blood pressure is fluctuating from beat to beat and is influenced by daily activities measurements can differ from measurement to measurement.

## Maintenance and storage

---




- For cleaning the monitor only use a soft, lightly moistened cloth.
- **Do not use petroleum spirits, thinners or similar solvents!**
- Stains on the cuff can be carefully removed using a moist cloth and soapy water.
- **Do not wash the cuff or make it wet. Please make sure no water will enter the main unit!**
- Keep your monitor in the storage case protected from extreme temperatures and humidity.
- Do not carry out repairs of any kind yourself. If a defect occurs consult your OMRON dealer or distributor.
- Please take the batteries out if the monitor is not used for a long period (app. 2 months or more).
- Electro Magnetic fields can influence the measurement. Do not use a portable phone near the unit.

Disposal of this product and used batteries should be carried out in accordance with the national regulations for the disposal of electronic products.

## Error indication

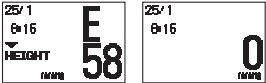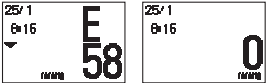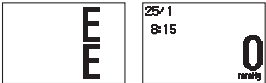If the measurement is not made correctly, the following error can be displayed.

Error	Cause	How to correct
 <p>After several seconds</p>	You moved your wrist during measurement.	Measure your blood pressure again without moving your wrist until the measurement is completed. (Refer to Page 10.)
 <p>After several seconds</p>	You moved your arm or body, or talked during measurement.	Measure your blood pressure again without moving your arm or body. (Refer to Page 10.)
 <p>After several seconds</p>	You moved your arm or body extensively during measurement.	Measure your blood pressure again without moving your arm or body. (Refer to Page 10.)
 <p>Battery replacement mark flashes.</p>	Batteries are worn out.	Replace both batteries with new ones. (Refer to Page 6.)
 <p>Only E and numerals are displayed.</p>	Trouble caused by abnormality in function.	Please contact the nearest OMRON dealer.

If you cannot resolve the problem even after taking the above-mentioned measures, please contact the nearest OMRON dealer.

# Troubleshooting

Trouble	Possible cause	How to correct
Nothing is displayed when you push the on/off button.	Batteries are worn out.	Replace the batteries with new ones. (See to Page 6.)
	The directions of polarities of inserted batteries (+ and -) are wrong.	Insert the batteries correctly. (Refer to Page 6.)
Measurement cannot be made. Or the blood pressure values are displayed abnormally high (low).	The wrist cuff is not held at the level of the heart.	Measure your blood pressure again in correct posture. (Refer to Page 9.)
	The wrist cuff is not wrapped firmly.	Wrap the wrist cuff correctly. (Refer to Page 7.)
	Too much stress is applied on your shoulder or arm.	Relax yourself and make measurement.
	You are talking or moving your hands during measurement.	Be still and make measurement.
The blood pressure values differ from the values measured at brachial.	Refer to the section on "Wrist Blood Pressure Monitor." (Refer to Page 18.)	
Blood pressure values are displayed before the monitor inflates to the required value.	Air is leaking from the wrist cuff.	Contact the store where the unit was sold or the nearest OMRON distributor.
Blood pressure value fluctuates each time of measurement. Or the blood pressure values are displayed abnormally high (low).	Blood pressure value tends to fluctuate depending on the mental state during measurement or the measuring time. Breathe deeply to relax yourself and make measurement again. (Refer to Page 19.)	

If you cannot resolve the problem even after taking the above-mentioned measures, please contact the nearest OMRON dealer.

## Spare parts and accessories

OMRON considers, due to the used high company and quality standards, the main unit as a non-serviceable part because of the necessity of proper calibration after replacement of high-tech components.

The cuff can be replaced by an authorised expert.

## Specifications

Product Description :	Automatic Blood Pressure Monitor
Model :	R7 (HEM-637-E2 / HEM-637-E7)
Display :	Dot matrix LCDdisplay
Measurement :	Oscillometric method
Measurement Range :	Pressure; 0 to 299 mmHg, Pulse rate; 40 to 180 /min.
Accuracy :	Pressure; within +/-3 mmHg, Pulse rate; within +/-5% of
reading	
Inflation :	Automatic inflation with pump
Deflation :	Automatic rapid deflation
Power supply :	Two LR03 alkaline batteries (AAA)
Battery Life :	Approximately 300 times of measurement when using alkaline batteries (In case of the monitor is used three times a day inflating to 170 mmHg at the room temperature of 22°C)

### Operating Temperature

and Humidity : +10° to +40°C, 30 to 85% RH

### Storage temperature

and humidity : -20° to +60°C, 10 to 95% RH

### Measurable

circumference of wrist : 135 to 215 mm

Weight of Main Unit : Approximately 150 g (not including batteries)

External Dimensions : 78 (W) x 65.5 (H) x 37.2 (D) mm (not including the wrist cuff)

Electric Shock Protection : Internal power supply appliance type B

Accessories : Two LR03 alkaline batteries for monitoring use, storage case.  
Instruction Manual, blood pressure pass, guarantee card.

• In the interest of product improvement, specifications are subject to change without notice.

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Made in China

This device fulfils the provisions of the EC directive 93/42/EEC (Medical Device Directive). This blood pressure monitor is designed according to the European standard EN1060, Non-invasive sphygmomanometers Part 1: General Requirements and Part 3: Supplementary requirements for electromechanical blood pressure measuring systems.

This OMRON product is produced under the strict quality system of OMRON HEALTHCARE Co., Ltd., Japan. The Core component for OMRON blood pressure monitors, which is the Pressure Sensor, is produced in Japan.

=Type B

C€0197

<b>Manufacturer</b> 	<b>OMRON HEALTHCARE Co., Ltd.</b> 24, Yamanouchi Yamanoshita-cho, Ukyo-ku Kyoto, 615-0084 JAPAN
<b>EU-representative</b> <div data-bbox="83 301 260 375" style="border: 1px solid black; padding: 5px; display: inline-block;"> <div style="border-right: 1px solid black; padding: 0 10px;">EC</div> <div style="padding: 0 10px;">REP</div> </div>	<b>OMRON HEALTHCARE EUROPE B.V.</b> Kruisweg 577, 2132 NA Hoofddorp THE NETHERLANDS <a href="http://www.omron-healthcare.com">www.omron-healthcare.com</a>
<b>Production facility</b>	<b>OMRON (DALIAN) CO., LTD.</b> Economic & Technical Development Zone Dalian 116600, CHINA
<b>Subsidiary</b>	<b>OMRON HEALTHCARE UK LIMITED</b> Opal Drive, Fox Milne Milton Keynes, MK15 0DG U.K.
	<b>OMRON MEDIZINTECHNIK HANDELSGESELLSCHAFT mbH</b> John-Deere-Str. 81a 68163 Mannheim, GERMANY <a href="http://www.omron-medizintechnik.de">www.omron-medizintechnik.de</a>
	<b>OMRON Santé France SAS</b> 14, rue de Lisbonne 93561 Rosny-sous-Bois Cedex, FRANCE

Made in China

<b>Manufacturer</b> 	<b>OMRON HEALTHCARE Co., Ltd.</b> 53, Kunotsubo, Terado-cho, Muko, Kyoto, 617-0002 JAPAN		
<b>EU-representative</b> <table border="1" data-bbox="76 265 252 339"> <tr> <td data-bbox="76 265 163 339">EC</td><td data-bbox="163 265 252 339">REP</td></tr> </table>	EC	REP	<b>OMRON HEALTHCARE EUROPE B.V.</b> Scorpius 33, 2132 LR Hoofddorp THE NETHERLANDS <a href="http://www.omron-healthcare.com">www.omron-healthcare.com</a>
EC	REP		
<b>Production facility</b> <b>Produktionsstätte</b> <b>Productiefaciliteit</b> <b>Site de production</b> <b>Oficina de producción</b> <b>Stabilimento di produzione</b>	<b>OMRON (DALIAN) CO., LTD.</b> Dalian, CHINA		
<b>Subsidiary</b> <b>Niederlassung</b> <b>Dochteronderneming</b> <b>Succursale</b> <b>Empresa filial</b> <b>Consociata</b>	<b>OMRON HEALTHCARE UK LTD.</b> Opal Drive, Fox Milne, Milton Keynes MK15 0DG, U.K.		
	<b>OMRON MEDIZINTECHNIK HANDELSGESELLSCHAFT mbH</b> John-Deere-Str. 81a, 68163 Mannheim GERMANY <a href="http://www.omron-medizintechnik.de">www.omron-medizintechnik.de</a>		
	<b>OMRON SANTÉ FRANCE SAS</b> 14, rue de Lisbonne 93561 Rosny-sous-Bois Cedex, FRANCE Uniquement pour le marché français: OMRON Service Après Vente N° Vert 0 800 91 43 14 <a href="mailto:consommateurs@omron-sante.fr">consommateurs@omron-sante.fr</a> <a href="http://www.omron-sante.fr/contact">www.omron-sante.fr/contact</a>		

Made in China  
Hergestellt in China  
Geproduceerd in China  
Fabriqué en Chine  
Fabricado en China  
Prodotto in Cina

The familiar way

Professional diagnostic instruments  
of the highest standards.

Product Overview

# Content

3	Preface
4-5	Riester today
6-11	Sphygmomanometers
12-17	E.N.T./Ophthalmic instruments Dermatological instruments
18	Vital signs Monitors Pulsoximeter
18	Thermometers
19-21	Stethoscopes
22-23	Head mirrors, Headlights, Examination lamps
24	Laryngoscopes
25	Gynaecological instruments
25	Percussion hammers
26	Tuning forks
26	Blood stasis products
26	Pulmonary pressure measuring devices
26	Dynamometers
27	Pressure infusion instruments
27	Veterinary instruments
27	Diagnostic Sets

## Diagnostic instruments and devices of the highest standards.

"Riester – The familiar way " stands above all for the trust that our customers have in the high quality and safety, the well-known precision and the pioneering innovation potential in almost every country of the world. Designed for demanding everyday use in doctors' practices and hospitals, many Riester products have been and still are being created through close cooperation between our experienced research and development department and practitioners from government research centres and university hospitals. This is, however, only one element of the Riester product philosophy, because in all our activities it is the human aspect which still remains the centre of our attention. Furthermore, we even guarantee our promise of quality and safety in writing. All of our products have the CE mark, a sign for special quality, and are supplied with the "Original Riester" certificate of authenticity – no product leaves our premises until it has passed exacting quality tests. "The familiar way" stands not only for customer confidence but also for the strength, flexibility and individuality of a long tradition company which is still being successfully managed. Conspicuous proof of this is given by the

ISO-9001 certificate which we were awarded on 30.10.1997 – the culmination of more than 50-years of company history in development and production of medical diagnostic instruments. A history of trust which we are proud of, but which also commits us in particular for the future. With our highly motivated specialist personnel, who is trained most extensively in house, we are able to successfully master the balancing act between our continuity and future-oriented modernity. And if you wish, this can also include you as a partner who places trust in the perfection and innovative power of traditional methods. We would be pleased to welcome you as a member and customer of the great Riester family, if this is not already done.

Yours very sincerely

Management and Staff  
Rudolf Riester GmbH

## Riester – success based on experience.

Since 1948 when Riester has been founded, the company and its employees have been committed to the engineering and production of diagnostic instruments of the highest standards.

Today Riester employs 150 people at its headquarters in Jungingen/Germany and distributes its comprehensive range of products in 150 countries.

The major product-lines are sphygmomanometers, E.N.T. and ophthalmic instruments, stethoscopes, laryngoscopes and examination lamps.

Quality – because it pays: Riester is today among the best reputed companies in the market of diagnostic instruments - worldwide.

## Always one step ahead – Riester instruments.

Not resting on one's laurels, that's one of the basic principles at Riester.  
Due to new product-developments of the previous years Riester performed growth rates above average:

→ 2003

**Stethoscope line –**  
verified acoustic  
performance.

→ 2004

**F.O. Laryngoscopes –**  
XL illumination performance and  
verified durability of F.O. blades.

→ 2005

**ri-former® -**  
the modular diagnostic station

→ 2006

**ri-focus® LED –**  
LED illuminated headlamp

→ 2007

**ri-scope® L –**  
Instruments for ENT and  
Ophthalmology with LED /  
Li-Ion technology

→ 2008

**ri-vital® spot-check monitor**  
precise, robust, time effective  
R1 **shock-proof**  
The new standard for aneroid blood  
pressure measurement


## Quality meets precision.

### Success of our products „Made in Germany“ doesn't happen by accident.

That's why the manufacturing of quality components and assembly of Riester sphygmomanometers is done in Germany. Thus resulting precision and product longevity have been proven by test results. Riester's latest inventions shock proof and stable needle deflection are going to be patented technologies.

Riester aneroid sphygmomanometers differ not only from precision comparable to mercury devices, but also from a nearly unlimited lifetime. Another proof of its unparalleled quality is the diaphragm pressure-loading capacity of up to 600 mm Hg. Choose from our large product variety with screw- or plug-in connections and a large variety of cuffs. Riester sphygmomanometers are also available as latex-free versions, with high precision air release valves and micro filters to protect the measuring system.

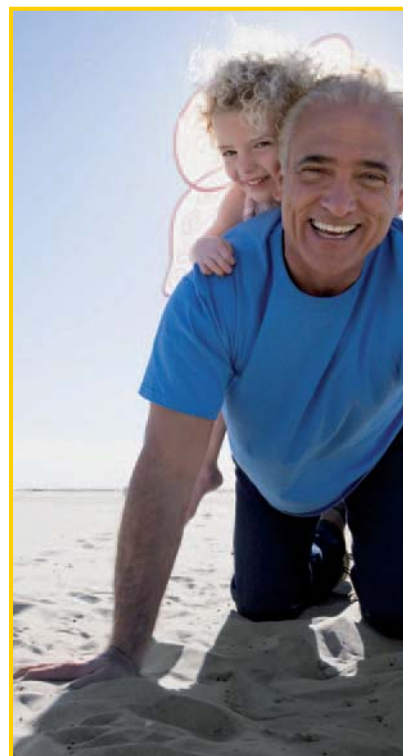SHOCK  
PROOF

### minimus® II

**Metal Aneroid sphygmomanometer. Single-tube model.**

- High-polished chromium-plated metal casing with metal ring, corrosion-protected.
- Metallic screw tube connection at the top of the manometer for undisturbed handling.
- Metal valve base.
- Spoon-shaped handle made of stainless steel.
- Linear aluminium scale Ø 49 mm, easy to read.
- Latex bulb with chromium-plated precision air-release valve, fine regulation and non-wearing.
- In vinyl zipper case.

### minimus® III

**Metal Aneroid sphygmomanometer. Twin-tube model.**

- High-polished chromium-plated metal casing with metal ring, corrosion-protected.
- Metallic screw tube connections.
- Metal valve base.
- Spoon-shaped handle made of stainless steel.
- Linear aluminium scale Ø 49 mm, easy to read.
- Latex bulb with chromium-plated precision air-release valve, fine regulation and non-wearing.
- In vinyl zipper case.

### R1 shock-proof

**R1 shock-proof. Exact precision that lasts a lifetime!**

- Patented (pending) shock-proof technology
- Patented (pending) air duct for stable pointer movements
- Fast, ergonomic inflation due to optimized ball design with integrated spoon
- Non-resistant opening of air release valve ensures precise and comfortable operation
- **5 years warranty on calibration**
- Shock-proof up to a height of 120cm. Conventional aneroid sphygmomanometers are shock-sensitive and its measurement accuracy may be impaired after an impact or a fall to the floor, so may no longer be used. This results in inaccurate readings or considerable logistic costs.
- Focus Green Silver 2008 Design Award
- Scale: Ø 52 mm white with blue pointer or black with red fluorescent pointer
- Valve base and air duct with metal precision valve
- Two-component ABS casing
- Metal threaded tube connection at the top of the casing for fast cuff exchange and ergonomic operation
- Maximum error tolerance of  $\pm 3$  mm Hg
- Extensive selection of various cuff sizes
- Incl. Vinyl zippered bag
- Weight: 173 gram
- Set with storage box and 3 cuffs (large adults, adults, children) and R1 shock-proof with white scale

### precisa® N

**Aneroid sphygmomanometer in aluminium- or ABS-plastic casing.**

- Metal valve base.
- Spoon-shaped handle made of stainless steel.
- Linear aluminium scale Ø 63 mm, easy to read.
- Latex bulb with chromium-plated precision air-release valve, fine regulation and non-wearing.
- In vinyl zipper case

#### Plastic casing:

- Sturdy ABS plastic casing.
- Metallic screw tube connection.

#### With aluminium casing:

- Sturdy and lightweight casing made of aluminium.
- Plug-in tube connection made of metal.

### babyphon®

**Suitable for children, all in blue. Palm-style sphygmomanometer for paediatrician's offices and hospitals.**

- Bear imprint on the cuff for children, infants and new-born children.
- Choice of manometer with chromium-plated metal casing or sturdy plastic casing.
- Ergonomic spoon-shaped handle made of stainless steel.
- Metal manometer with aluminium scale Ø 49 mm, plastic manometer with aluminium scale Ø 63 mm, both linear and easy to read.
- Latex bulb with chromium-plated precision air-release valve, fine regulation and non-wearing.
- In padded nylon zipper case.
- Set II offers a combination including the paediatric stethoscope duplex® baby in a matching colour.

## exacta®




Classic aneroid sphygmomanometer.  
Two-tube model.

- Choice of high-polished chromium-plated metal casing, black lacquered or black ABS-plastic casing.
- Stainless metal clamp for hooking the manometer to the cuff.
- Latex bulb with precision air-release valve, fine regulation and non-wearing.
- Linear aluminium scale Ø 49 mm, easy to read.
- In vinyl zipper case.

## ri-san®

Palm-style sphygmomanometer in consistent ergonomic design.

- Regulation of the air-release valve by simple finger-tip control.
- Plug-in tube connection at the top of the manometer for undisturbed handling.
- Latex-free bulb for inflating the cuff.
- Linear aluminium scale Ø 64 mm, easy to read.
- In vinyl zipper case.
- Model colours: slate grey, blue, saffron and green.
- Suitable for left and right-handed persons.
- Sturdy ABS-plastic casing.
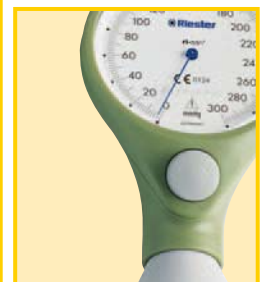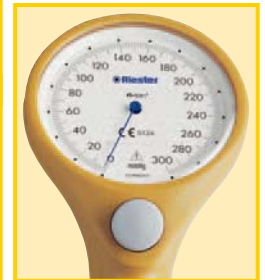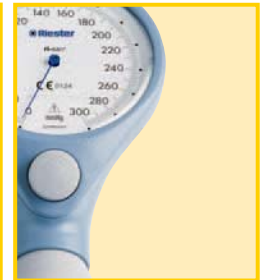
## ri-san® combi set

The special Riester service:  
Practical combi-sets with unique attractive combined price advantage!

- Ideal for physicians and students and for every hospital.
- The sets contain instruments in matching colours.
- colours: slate grey, blue, saffron, green.
- Individual choice of product combination on demand.

## sphygmotensiophone

## fortelux®H

## pen-scope®

## ri-mini®

## duplex®

Aneroid two-tube sphygmomanometer in metal casing.  
The chromium-plated pressure gage housing fits firmly on the cuff.

- High-quality cotton cuff finished in artificial leather with a metal clip to pull through the cuff. Sizes available for adults, children and infants.
- Latex bulb with precision air-release valve, fine regulation and non-wearing.
- Linear aluminium scale Ø 49 mm, easy to read.
- In padded nylon zipper case.

## big ben®

### Exact precision that lasts a lifetime.

- Available as desk, wall, floor or anaesthesia models, either round or square.
- Generously dimensioned scale (Ø 147.2 mm or 145.5 mm) and high-contrast scale layout for optimum readability.
- Spacious basket on the reverse side for tubes and cuff.
- Latex bulb to inflate the cuff.
- Spiral tube can be extended to 3 m.
- Non-wearing air-release valve with fine adjustment.
- Precision mechanism attached to a special support to ensure complete stability of the measuring system and the zero-point.

#### Wall model:

- Ideal for every hospital room.
- Can be swivelled horizontally 65° to the left and right.
- Glass-fibre reinforced wall bracket for maximum strength (comparable to steel).
- Simple conversion to the rail model.

#### Anaesthesia model:

- Same features as the wall model.
- Simple conversion to the rail model.

#### Floor model:

- Choice of two firm stands:
  - Ø 38 cm, powder-coated steel, five casters or
  - Ø 55 cm, antistatic, stainless steel, five casters, two locking casters for completely safe standing.
- Stand height of 1.05 m adjustable to 1.30 m.

Acquiring mercury sphygmomanometers is a matter of confidence. The Riester quality-manufacturing is making the difference. Riester sphygmomanometers differ from others not only by highest precision and safety, but also by nearly unlimited lifetime – since 1948. All devices are equipped with micro filters to protect the measuring system, special seals against mercury contamination and mercury locks.

### diplomat-presameter®

**High-precision mercury sphygmomanometer. Desk model in a metal case.**

- Precision glass column with internal diameter of 4.2 mm  $\pm$  0.2 mm.
- Chromium-plated precision air-release valve, non-wearing and with fine regulation.
- Nickel-plated metal parts.
- Latex bulb with air-release valve.
- Cuff with two-tube latex bladder.
- Unbreakable, long-life metal case with large compartment for cuff and latex bulb.
- Model colours: light grey, slate grey, blue, saffron and green.

### nova-presameter®

**High-end mercury sphygmomanometer. Desk, wall or floor model in metal casing.**

- Precision glass column with internal diameter of 4.2 mm  $\pm$  0.2 mm. The floor model has an internal diameter of 5.0 mm and tilted scale for the best legibility.
- Chromium-plated precision air-release valve, non-wearing and with fine regulation.
- Chromium-plated metal parts.
- Desk model including bulb holder.
- Latex bulb with air-release valve.
- Cuff with two-tube latex bladder.
- Spiral tube of floor and wall models can be extended to 3 m.
- Unbreakable, durable metal case with two-coloured lacquering light grey/black with large compartment for cuff and latex bulb.

### nova-presameter®

**Mobile mercury sphygmomanometer.**

- Precision glass column with internal diameter of 5.0 mm  $\pm$  0.2 mm and tilted scale for the best legibility.
- Chromium-plated precision air-release valve, non-wearing and with fine regulation.
- Chromium-plated metal parts.
- Latex bulb with air-release valve.
- Cuff with two-tube latex bladder.
- Spiral tube can be extended to 3 m.
- Unbreakable, durable metal case with two-coloured lacquering light grey/black with large compartment for cuff and latex bulb.
- Choice of two firm stands:
  - Ø 38 cm, powder-coated steel, five casters or
  - Ø 55 cm, stainless steel, five casters, two locking casters for completely safe standing.
- Stand height of 1.05 m adjustable to 1.30 m.

### empire®N

**Mercury sphygmomanometer of ABS-plastic casing.**

- Attractive design combined with sturdy everyday-use quality, available as desk, wall, floor or anaesthesia model.
- Precision glass column with internal diameter of 4.2mm,  $\pm$  0.2mm.
- Chromium-plated precision air-release valve, non-wearing and with fine regulation.
- Chromium-plated metal parts.
- Latex bulb with air release valve.
- Cuff with two-tube latex bladder.
- Spiral tube can be extended up to 3m.

### global de luxe®

**Wall mounted mercury sphygmomanometer with swivelling metal casing.**

- Metal casing can be swivelled 50° to the left and right.
- Precision glass column with internal diameter of 4.2mm  $\pm$  0.2mm
- Chromium-plated precision air-release valve, non-wearing and with fine regulation.
- Chromium-plated metal parts.
- Latex bulb with air-release valve.
- Cuff with two-tube latex bladder.
- Spiral tube can be extended up to 3m.

### sanaphon®

Classic palm-style sphygmomanometers for precise self-monitoring. Stethoscope included.

- Suitable for individuals with heart rhythm disorders.
- Stethoscope chest-piece permanently attached to the cuff.
- Aluminium casing with metal ring, protected against corrosion.
- Plug-in tube connection made of metal.
- Metal valve base.
- Spoon-shaped handle made of stainless steel.
- Linear aluminium scale Ø 64 mm, easy to read.
- Latex bulb with chromium-plated air-release valve, fine regulation and non-wearing.
- Also available with precisa® N plastic/metal or ri-san® manometers.
- In padded nylon zipper case.
- Personal blood pressure record booklet included.
- In attractive promotional cardboard display box.

### ri-san®

Palm-style aneroid sphygmomanometer for exact self-monitoring with stethoscope.

- Ergonomic design in attractive colours.
- Suitable for right and left-handed individuals.
- Suitable for persons with heart rhythm disorders.
- Stethoscope chest-piece permanently attached to the cuff.
- Air-release by simple finger-tip control, fine regulation.
- Plug-in tube connection.
- Latex-free bulb.
- Linear aluminium scale Ø 64 mm, easy to read.
- In padded nylon zipper case.
- Personal blood pressure record booklet included.
- In attractive promotional cardboard display box.

### ri-champion® N

ri-champion® N: Fully automatic digital blood pressure monitor for precise and fast measurement.

The clinically validated device according to BHSstandard A/A offers beside an easy handling and irregular heartbeat recognition also 99 memories, pulse tone, date and time. It operates with AA-batteries (included) or an optional mains adapter.

- Measurement of blood pressure and pulse
- Meets the high BHS-standard A/A
- Recognition of arrhythmia
- Large display
- Optical and acoustical pulse signal
- 99 memories
- Date and time
- Cuff storage
- Optional mains adapter

## ri-scope® L – The new era of light and optics.

The ri-scope® L opens a new era of high-performance and user-friendly diagnostic instruments. Li-ion/LED technology and new high-performance instruments significantly reduce costs while improving both effectiveness and reliability. Instruments for E.N.T. and Ophthalmology.


### Lithium-ion battery – The power plus.

Memory effect and self discharge are not match: the light-weight lithium-ion battery used in the ri-scope® L has an unsurpassed lifetime – with higher loading capacity and sustained high-power performance, even after well over 1000 charge cycles.

### LED technology – performance for life.

Nearly unlimited lifetime with the highest luminosity: these advantages make the LED technology used with the ri-scope® L so unique. Moreover, the lamp inside the ri-scope® L never gets hot, provides true-colour, even illumination of the area under examination and is practically guaranteed never to fail, in contrast to conventional halogen or xenon lamps.
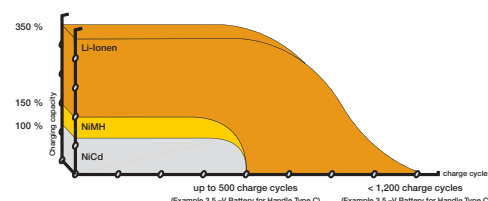
### rheotronic® – Patented technology for LEDs and bulbs.

Handy Riester handles with continuously adjustable **rheotronic®**-for regulating light intensity are the basis for a safe examination. With automatic security switch-off.

### ri-scope® L Otoscopes/ Ophthalmoscopes – New standards for eye and ear diagnostics!

These completely new instruments were developed with input from doctors to set new standards for optical quality, ergonomics and performance options.

**LED + Li-Ion**  
*inside*

Otoscope with battery handle Type AA

Ophthalmoscope with battery handle Type C


## First-class ENT and ophthalmic instruments.

### ri-scope® L.

ri-scope® L with LED, Xenon or Halogen illumination is the first choice for the challenges in E.N.T and ophthalmology – everyday. Its diverse combination possibilities and technical options make it a highly valued instrument in every doctor's practice and clinic. All ri-scope® L instruments heads are equipped with a bayonet fitting for quick and secure attachment to the handle. A large selection of energy sources and practical charger units is available.

ri-scope® L is developed and manufactured in Germany.

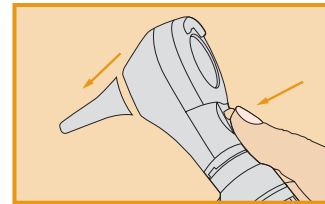The SpecEjec® device enables touchless ejection of used ear specula for hygienic and simplified handling.

### ri-scope® L Otoscop

#### ri-scope® L 1, 2 and 3

- ri-scope® L Otoscope with optical glass and reduced reflection.
- Basic model L1 with reflection minimized direct illumination Halogen HL 2.5V or Xenon XL 3.5V.
- Bi-directional swivelling optical glass with 3-fold magnification.
- Suitable for pneumatic tests (supplied without ball).
- Available with optional C or AA handles.
- L2 and L3 version with fiberoptic for optimal beaming and transmission of the light.
- High-end version L3 equipped with operation lens with 4-fold magnification and specula ejection device SpecEjec®.

### ri-scope® L Ophthalmoscop

#### Corrective lenses L1/L2

D+	1	2	3	4	5	6	7	8	9	10	12	15	20	40	
D-	1	2	3	4	5	6	7	8	9	10	15	20	25	30	35

Corrective lenses from +40 through -35 diopters guarantee a sharp picture, even in cases of especially severe refraction anomalies.

#### ri-scope® L 1, 2 and 3


- All ri-scope® L ophthalmoscopes with modular design impress with their high-performance optics with aspherical condenser lens and reduced reflection.
- The new ri-scope® L ophthalmoscopes have a thumbrest for securing the ophthalmoscope, which enables easy control of all the elements with one finger.
- Reflexes from the cornea and iris are avoided by separating the observational and illuminating light beams (Gullstrand's principle).
- 3 different versions L1, L2 and L3 available with correction lenses with up to 89 diopter values
- Easy to operate aperture hand-wheel for apertures and filters included in all three versions.
- High-end version L3 includes a focussing device with zero retention for force for the quick correction of ametropias.

### ri-scope® L sets


Wide range of ri-scope® sets with a cost-effectiveness which is particularly convincing. Batteries and charging units have to be ordered separately. Individual combinations of instruments on demand.

- 2.5V instruments including HL halogen lamps
- 3.5V instruments including XL xenon or LED lamps
- All sets are supplied in a hard case with safety catches including spare lamps
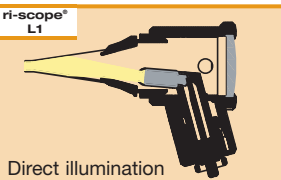
#### Direct illumination

Reflection reduced, cost effective direct illumination with 2.5 V HL halogen or 3.5V XL xenon lamps.
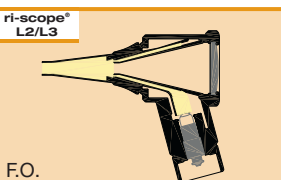
#### F.O.

Fiber optics with 2.5V HL halogen or 3.5V LED illumination ensure an unhindered, glare-free view.

#### Apertures L1

#### Apertures L2

#### Filters

#### Apertures L3

#### Filters

## ri-former® – the modular diagnostic station.

### ri-former® is adaptable to your needs. Each time anew.

The modular diagnostic station ri-former® with LED illumination can be customized according to your needs, it can be expanded to up to 5 handles. ri-former® always perfectly matches the requirements in clinic and medical practice.

Available with or without optional digital clock. The constant output voltage independent of the power supply prolongs the service life of the lamps and guarantees optimal light intensity.

The handles automatically switch on or off when removing or replacing them in the handle unit, while the adjustable light intensity setting is maintained until the next examination.

Instrument heads are available with or without anti-theft device.

Up to 5 handle modules plus ri-spec® specula dispenser plus the big ben® square sphygmomanometer, as well as the ri-thermo® N wall unit can be mounted on the right or left side of the ri-former® basic module.

**LED**  
*inside*

## pen-scope®

### Standard Pocket Oto-/Ophthalmoscopes

- Choice of otoscopes, ophthalmoscopes or oto-/ophthalmoscopes in attractive nylon case.
- Highly attractive cost-effectiveness.
- Practical metal clip on the handle for attaching to the doctor's coat.
- Otoscope with direct illumination, alternatively Standard- or Halogen lamps, 2.5x magnifying glass, removal for inserting external instruments.
- Ophthalmoscope with four different apertures and correction lenses from 0 to +20 and 0 to -20 diopters, easy to operate focusing wheel and integrated eyeglass protector.

D+	1	2	3	4	6	10	15	20
D-	1	2	3	4	6	10	15	20

### pen-scope® Otoscope

- 2.7V vacuum lamp or 2.5 HL halogen lamp.
- Direct illumination, with an upgraded, reflection-minimized illumination.
- Removable magnifying lens with 2.5-fold magnification.
- Possibility to insert external instruments when magnifying lens is removed.
- Screw fitting for secure attachment to the handle.
- Sturdy adjustment ring to turn the instrument head to the ideal position to the handle.
- Sealed system for performing pneumatic otoscopy.
- Extremely impact resistant casing made of glass-fibre reinforced plastic.
- Simple exchange of lamp at the base of the instrument head.
- Including ear specula 4mm and 2.5mm.

Large circle

Small circle

Semi circle

Fixation star


## ri-mini®

D+	1	2	3	4	6	10	15	20
D-	1	2	3	4	6	10	15	20

### Halogen/F.O. illuminated Oto-/Ophthalmoscopes

- 2.5V F.O. and Halogen illumination.
- Selection of otoscopes, ophthalmoscopes or oto-/ophthalmoscopes in hard case with safety catches.
- Highest cost effectiveness.
- Metal clip on the handle for attaching to the physician's coat.
- Otoscope with fibre-optics system for an optimum light beam and thus the best possible lighting results.
- Magnifying lens with 2.5-fold magnification, removable, ideal for introducing external instruments.
- Ophthalmoscope with four different apertures and correction lenses from 0 to +20 and 0 to -20 diopters, easy to operate focusing wheel and integrated eyeglass protector.

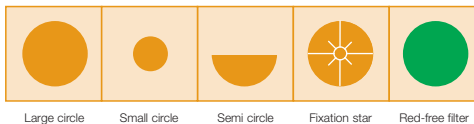Large circle

Small circle

Semi circle

Fixation star

Red-free filter

## uni®

### uni® Oto-/Ophthalmoscope.

#### uni® I Otoscope:

- Otoscope head, chromium-plated, sturdy.
- With plug-in connector.
- Removable cover glass.
- Chromium-plated removable lens with 4x magnification.
- With Vacuum lamp 2.7V , 2.5V HL Halogen or 3.5V Xenon lamps.

#### uni® II May-Ophthalmoscope:

- Correction lenses from 0 to +20 and 0 to -20 dioptries.
- Aperture: large circle.
- With Vacuum lamp 2.7V , 2.5V HL Halogen or 3.5V Xenon lamps.

#### uni® III Otoloscope/May-Ophthalmoscope:

- All models are fitted with battery handle type C for two batteries type C, and with rheostat for regulating the light intensity.

## econom®

### The complete diagnostic set.

- Vacuum lamp 2.7V , 2.5V HL Halogen or 3.5V Xenon lamps.
- A choice of two sets is available: with expandable nasal speculum or with Ø 9 mm nasal speculum.
- Alternatively supplied with halogen or xenon lamps and rechargeable handles with NiMH batteries.
- Apart from uni® otoscope and May-ophthalmoscope included are: Nasal speculum, bent arm illuminator, 2 laryngeal mirrors, tongue blade holder, 3 ear-specula, handle and spare lamp.
- Nasal speculum expandable, chromium-plated metal.
- Bent arm illuminator: Chromium-plated metal, for standard wooden tongue blade sizes or laryngeal mirrors , with battery handle type C for two batteries type C.
- In sturdy black case.

## fortelux® H

### High performance diagnostic penlight. The ideal addition to ri-mini® and pen-scope® instruments!

- Bright concentrated halogen light 2.5V.
- Handle AA for two batteries type AA 1.5V.
- Practical metal clip for attaching to the physician's coat.
- Including tongue blade holder for commercially available wooden tongue depressors.
- In modern nylon case with belt loop.
- Attractive choice of colours, matching ri-mini®, pen-scope® as well as other Riester products.

## fortelux® N

### Economically priced diagnostic penlights. Metal casing in attractive colours.

- With 2.2 V vacuum lamp.
- Including 2 batteries type AAA.
- Aluminium outer casing with disinfectable special lacquer.
- With practical metal clip for attaching to the physician's coat.
- Simple on/off function through contact with the metal clip.
- Supplied singly or in packs of 6 (assorted colours or six per colour).

## ri-light®

### Cost-effectively priced diagnostic set for students.

- With diagnostic penlight fortelux® N, two batteries type AAA included.
- Tongue blade holder made of plastic.
- Laryngeal mirrors made of plastic.
- Blue filter in white plastic mounting.
- In plastic case.

## ri-vision®

### ri-vision® with slit-light retinoscope or spot-light retinoscope and L2 ophthalmoscope.

- Available in 2.5V HL halogen or 3.5V XL xenon illumination.
- In slit- or spot- version available.
- Two fixation cards for dynamic retinoscopy.
- 2 spare lamps included.

## ri-scope® Retinoskop

For recognition of myopia and hyperopia as well as astigmatism.

- XL 3,5V xenon lamp or HL 2,5V halogen lamp.
- Two models are available: The spot-light retinoscope projects a circular light beam. The slit-light retinoscope with a light beam in the form of a line, simplifies recognition and the determination of astigmatic refractive errors.
- Simple operation with knurled thumb screw. The line and spot image can be focussed with the operating element and turned 360°, angle can be read off the integrated scale.
- Holder for hanging and fixing the fixation cards into position for dynamic retinoscopy.
- Two fixation cards and one spare lamp are supplied within the set.
- Integrated eyeglass protection.
- Bayonet fitting for fast and secure attachment to the handle.
- Dust-tight, very sturdy and light casing made of impact resistant plastic.
- Simple exchange of the lamp at the base of the instrument head.

## schiotz tonometer

Highly precise eye tonometer according to Prof. Schiötz for measuring intra-ocular pressure.

- All vital parts are made of stainless steel.
- High quality agate bearing for an extremely long service life.
- Perfect reading of the scale with red pointer.
- Supplied with 3 weights (5.5g, 7.5g, 10g) and a conversion table.
- Complete in a black hard case with velvet-look inserts.
- Models with straight and inclined scale available.

## ri-derma® LED

Direct light microscopy for examining the skin.

- Particularly suitable for the differential diagnosis of pigmented skin changes, such as e.g. malignant melanoma.
- The upper skin layers are made visible by wetting with immersion oil. Surface structures, colour, pigment distribution, non-melanocytic structures can be assessed.
- Choice of HL® halogen lamps 2.5V and 3.5V LED lamps.
- Instruments (handles and heads) made of chromium-plated metal, sturdy and durable.
- Instrument head constructed to be dust-proof.
- Focusable magnifying glass with 10-fold magnification for the most precise diagnoses and detailed images.
- Two attachments which do not harm the skin – one with a scale from 0 to 10 mm for measuring skin changes and one without a scale; both can be disinfected.
- Eyeglass protector.
- Bayonet fitting for fast attachment of the instrument head to the handle.
- Simple exchange of the lamp at the base of the head.
- All handles and charging options from the ri-scope® are available.
- Comprehensive operating instructions.

## All Pupil II

WIRELESS

Unique performance and operating comfort

**High-contrast, clear images with the world's leading Keeler optics.**  
Riester's All Pupil II uses specially coated and individually laser-cut lenses for best sight – every day.

**Light and comfortable**

Riester's indirect ophthalmoscope only weighs 520 gram and is equipped with a particularly soft, cushioned and comfortable headband.

**Simple operating elements**

All Pupil II only has a unique single control to change aperture size and position in one movement to view through all pupil sizes.

**IR and UV filter**

An IR and UV filter is integrated for the patient's best possible protection.

**Easy and fast bulb replacement**

Easily accessible xenon bulbs enable rapid bulb replacement for uninterrupted work.

Cable connected or wireless power supply (patent pending)

## Binocular loupes

NEW!

Designed for comfort and performance

An enlarged view simplifies your work. What you can't see you can't work on. The advantages of Riester's binocular loupes are obvious: Improved working posture, reduced strain on the eyes and an even clearer view with K-LED lighting. Ideal for ocular and dental examinations as well as plastic surgery and veterinary professionals.

**Patient comfort**

More distance during treatment makes patients feel visibly better and more relaxed.

**Several magnification grades**

Riester magnifiers perform in accordance with your requirements. Multiple magnification levels enable you to see more, spot problems earlier and offer even higher levels of precision treatment.

**Attachment options**

Select between the light Riester spectacle frame or the brow band system for the use of your magnifiers.

ri-vital®

NEW!

ri-vital® spot-check monitor, engineered for the mobile, precise and fast recording of vital signs data at clinics, hospitals and outpatient treatments. Ergonomical design and highly developed technology resulting in reliable measurement results and time efficiency.

- NIBP, Pulse, MAP, SpO2, Temperature (optional)
- Special mode for neonates (blood pressure)
- Infrared ear thermometer with probe sheath dispenser (optional)
- Large, bright LCD display
- Simultaneous measurement of different vital parameters within seconds
- Shock resistant
- Easy and fast operation due to the simple menus and perfectly integrated storage options for sensors and cuffs.
- Efficient Li-Ion battery pack
- Portable and light
- Mobile stand (optional)
- 2 years warranty

ri-fox

Mobile precision in the smallest possible format! Monitor and sensor in a single device.

- Simple and fast measurement by attaching to the finger.
- Fast and precise measurement of the pulse and the oxygen content of the blood (%SpO2).
- The measurements are easy to read on the LCD display with optional background illumination.
- With acoustic pulse indicator.
- LCD-display of %SpO2, perfusion and pulse.
- Acoustic pulse indication.
- Battery status indication.
- Low battery consumption (approx. 24,000 measurements).
- commercially available batteries included, easy to change.
- Automatic on/off switch.

ri-thermo® N

Multifunctional Infra-red thermometer for professionals and private consumers – with CE 0124 approval, clinical validation.

- Three measurement functions: Apart from body temperature, exact and simple measurement of the temperature of liquids (e.g. baby bath-water or baby food) and the air temperature.
- Measurement takes approx. 1-2 seconds.
- Digital display for optimum readability of the measurement values in °C or °F.
- Memory function for 12 measurements.
- Probe covers: A new probe cover should be used for every measurement in the ear for reasons of hygiene and accuracy. No ear-measurement should be made without a probe cover!
- Two alkaline manganese batteries type AA 1.5V included for approx. 1,500 measurements.
- Supplied with 25 probe covers in attractive sales packaging.
- Probe covers are available in packs of 25, 100 or 1,000.

ri-gital®

Easy to use, robust and precise digital thermometer.

- For oral, rectal and under-arm measurement of body temperature.
- Highly sensitive temperature sensor for fast, precise measurements.
- Three-digit display.
- Memory function.
- Signal tone.
- Unbreakable and 100% water-tight.
- Measurement range: 32 °C to 43,9°C.
- Operating temperature: 10 °C to 40 °C (50 °F – 104°F).
- Accuracy: ± 0.1 °C.
- Automatic switch-off after 10 minutes.
- Includes long-life button cell battery for approx. 200 hours.
- In plastic storage box.

## Riester stethoscopes – There's simply more to hear.

### Outstanding acoustic performance in all frequency ranges.

Riester stethoscopes ensure perfect auscultation. „Made in Germany“ stands for selected materials and first-class workmanship to ensure perfect auscultation and outstanding durability. Current test protocols show: The newest Riester generation of stethoscopes has unique acoustic characteristics in comparison to all other models measured. Extra-soft, replaceable ear-tips with metal screw thread ensure a comfortable and effective sealing of the external auditory canal. Particularly flat and soft non-chill rims for improved adaptation on the skin are also included.

All Riester stethoscopes are supplied with replacement parts and comprehensive instructions for use.

24 months guarantee on all Riester stethoscopes.

### cardiophon: frequency analysis

20 Hz – 100 Hz max. 5,9 dB at 100 Hz



100 Hz – 1.500 Hz max. 28,3 dB at 320 Hz

### cardiophon

#### cardiology-stethoscope cardiophon. The Riester masterpiece.

- Especially in the frequency range between 200Hz and 500Hz, which is particularly important for the auscultation of heart sounds, the cardiophon provides unique acoustic performance.
- The cardiophon is suitable for the precise auscultation of adults and children.
- Double-headed chest-piece made of stainless steel.
- Special membranes on both sides for precise acoustics, Ø 48 mm and 36 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Anatomically shaped binaural made of stainless steel with reinforced, internal multiple spring.
- Soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Extra-strong Y-tube made of PVC with separate sound-conducting pathways for the left and right ears with particularly good audio-conductivity.
- Overall length: 72 cm.
- Latex-free.
- The cardiophon is supplied with two pairs of replacement ear-tips, a replacement membrane each and a name plate.

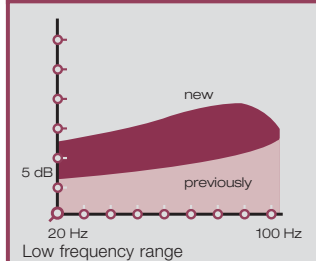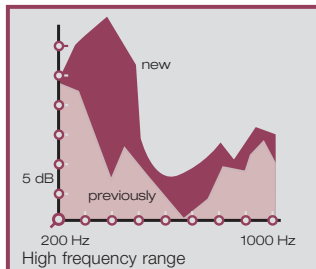Riester special membranes are extra-ordinary durable and permit perfect auscultation.

Vital signs Monitors  
Pulsoximeter

Thermometers

Stethoscopes

## duplex® de luxe

**Stethoscope in a class of its own for the auscultation of adults and children.**

- Considerably improved acoustics in all relevant frequency ranges.
- Double chest-piece made of either stainless steel or pleasantly light aluminium.
- Special membranes on both sides, Ø 48 mm and 36 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Matt chromium-plated binaural with integrated reinforced double spring.
- Extra-soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Extra-strong Y-tube made of PVC for shielding against undesired stray sounds.
- Latex-free.
- Overall length: 77 cm.
- duplex® de luxe is supplied with two pairs of replacement ear-tips, a replacement membrane each and a name plate.

## duplex®

**Quality stethoscope with double chest-piece. Best quality at an attractive price.**

- Considerably improved acoustics in all relevant frequency ranges.
- Double-headed chest-piece available in two types: chromium-plated brass or specially lightweight aluminium, galvanized.
- Special membrane for perfect auscultation, Ø 48 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Bell, Ø 36 mm.
- Chromium-plated binaural, which can be adjusted by turning. With external spring.
- Soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Strong Y-tube made of PVC to shield against undesired stray sounds.
- Latex-free.
- Overall length: 77 cm.
- duplex® is supplied with a pair of replacement ear-tips and a replacement membrane.

## tristar®

**The high quality stethoscope tristar® with three different sized double chest-pieces for adults, children/babies and new-born children is outstandingly suitable for perfect auscultation in all frequency ranges.**

- Considerably improved acoustic performance in all areas of application.
- Double chest pieces for duplex®, duplex® baby and duplex® neonatal made of aluminium with rapid-exchange connector.
- Special membranes, Ø 48 mm, 36 mm, 28 mm, for first-class acoustics.
- Particularly flat and soft non-chill rims for improved adaptation on the skin.
- Bells Ø 36 mm, 28 mm, 24 mm.
- Matt chromium-plated binaural with integrated reinforced double-spring.
- Optimal sealing of the external auditory canal with soft, replaceable ear-tips with metal screw thread.
- Particularly strong Y-tube made of PVC for optimal acoustics to shield against undesired stray sounds.
- Overall length: 79 cm, 78.5 cm, 78 cm.
- Latex-free.
- tristar® is supplied with two pairs of replacement ear-tips, a replacement membrane each and a name plate.

## anestophon®

**Innovative universal stethoscope with ergonomic flat chest-piece for particular comfort when measuring blood pressure. High quality standard at a reasonable price.**

- Considerably improved acoustics.
- Ergonomic handling thanks to novel design: extra-ordinarily flat aluminium chest-piece without edges, easy to push under the cuff when measuring blood pressure.
- Special membrane, Ø 48 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Matt chromium-plated binaural, which can be adjusted by turning.
- Soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Reinforced Y-tube made of PVC to shield against disturbing stray sounds.
- Overall length: 77 cm.
- Latex-free.
- The anestophon® is supplied in a cardboard display box with a pair of replacement ear-tips and a replacement membrane.

### duplex® baby

**Quality stethoscope with double chest-piece for children and babies.**

- Considerably improved acoustics.
- Particularly light precision chest-piece made of aluminium.
- Special membrane, Ø 36 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Bell, Ø 28 mm.
- Matt chromium-plated binaural, which can be adjusted by turning. With external spring.
- Soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Strong Y-tube made of PVC to shield against undesired stray sounds.
- Latex-free.
- Overall length: 76.5 cm.
- duplex® baby is supplied with a pair of replacement ear-tips and a replacement membrane.

### duplex® neonatal

**Quality stethoscope with double chest-piece for new-born children and infants.**

- Considerably improved acoustics.
- Particularly small precision chest-piece made of aluminium especially for new-born children and infants.
- Special membrane, Ø 28 mm, with particularly flat and soft non-chill rims for improved adaptation on the skin.
- Bell, Ø 22 mm.
- Matt chromium-plated binaural, which can be adjusted by turning. With external spring.
- Soft, replaceable ear-tips with metal screw thread for comfortable and effective sealing of the external auditory canal.
- Strong Y-tube made of PVC to shield against undesired stray sounds.
- Latex-free.
- Overall length: 76 cm.
- duplex® neonatal is supplied with a pair of replacement ear-tips and a replacement membrane.

### ri-rap®

**The universal stethoscope ri-rap® with double chest-piece and double tube system convinces by its flexible range of applications.**

- Very good acoustics in all frequency ranges.
- Strong double tube with separate sound-conducting pathways for the left and right ears guarantees extraordinary audio-conductivity.
- Alternatively with 40 cm or 80 cm tube length.
- Heavy chromium-plated double chest-piece with two membrane sides, Ø 45 mm and 32 mm.
- Special membranes with surface structure to adapt better to the skin.
- Both membrane sides can be simply and quickly refitted with the enclosed bells (Ø 33.5 mm, 27.3 mm, 20.3 mm) with the screw connections.
- Replacement parts box: three bells, each with a transparent membrane, two pairs of replacement ear-tips.
- Latex-free.
- ri-rap® supplied complete with replacement parts in an attractive cardboard display box with foam insert.

### pinard

**The pinard stethoscope is used to auscultate the heart sounds of the embryo in the womb.**

- Available in two models: sturdy plastic or particularly lightweight, polished aluminium.
- Supplied in a cardboard box.

Bright light for  
safe diagnosis.

LED  
inside

**Riester head mirrors and headlights show their particular quality in everyday use.** Head mirror and headlights have to provide optimal illumination, but also have to fit properly. The high level of wearing comfort means non-fatiguing operation even after long periods of use. Innovative light techniques guarantee optimum viewing conditions. Riester head mirrors, headlights and examination lights are developed and manufactured in Germany.

### ri-focus LED

#### Powerful head light with focusable LED lamp.

- Headband version, with built-in battery compartment and integral charge jack, can be optionally fitted with 2 lithium batteries size CR 123 A, 4 batteries size AAA (alkaline) or 4 rechargeable batteries size AAA (NiMH).
- Headband cleaning made easy by inner, removable and washable padding.
- Particularly comfortable with balanced, infinitely adjustable headband.
- Use without any hindering cables.
- On/off switch on battery compartment.
- Charge jack for plug-in charger integrated in battery compartment.
- Focusable LED lamp at a distance of 40 cm approx. Ø 40 mm - Ø 200 mm, infinitely adjustable.
- White LED, 1 Watt, 4.200° Kelvin, 6V.
- 7.100 Lux at 500 mm, 17.500 Lux at 300 mm.
- Battery operation: approx. 30 hrs. with lithium batteries, approx. 90 min. with alkaline batteries Type AAA.
- Rechargeable battery operation: approx. 90 min. with fully charged battery.
- Weight: 210 g without batteries, 242 g with CR 123 A lithium batteries, 259 g with rechargeable AAA batteries.
- Delivery in carton.

### ri-focus/ ri-focus reverse (halogen)

#### Optimum non-glaring illumination for all examinations.

- Halogen lamp 8 V for white and even light.
- The lens at the front of the light head can be focused to allow exact adjustment of the size of the illuminated examination area.
- Integrated heat sink for patient-friendly use.
- Step switch on cable for regulating light intensity in two steps.
- Electrical cable runs along the outside of the headband to permit unhindered use.
- Fast exchange of light heads using special fixture located on the side of the forehead.
- Ball joint for precise adjustment of the lamp head.
- Simple adjustment and fixation of the headband with help of the rotating aluminum knob.
- Soft, exchangeable foam protection around the forehead for comfort in wearing the headband.
- Rapid plug-in connector for attaching the headlight cable to the transformer ensures fast implementation.
- Available with 230 V or 120 V transformers.
- 7.000 Lux at 500 mm, 20.500 Lux at 300 mm
- Supplied in a sturdy plastic case.

### ri-focus CL (coldlight)

#### Powerful surgical head light.

- Precision optics with micro-focus and ball joint for a homogenous and precisely adjustable light point according to size and direction.
- Glass fibre cable for cold-light supply with centred cable duct for undisturbed working.
- Ergonomic headband with head bow and removable leather padding for comfortable use even over several hours, infinite height and width adjustment with central wheels adjustable with one hand.
- Optimum light output by choice of 3 different projectors:
  - HL 150 (150 W halogen)
  - XL 24 (24 W xenon)
  - XL 100 (100 W xenon).
- Fast push connector for connecting glass fibre cable and light source.
- Adapter for projectors of other quality manufacturers.
- Light exit: focusable, double lens system, ball joint.
- Glass fibre cable: length 220 cm, diameter 3.5 mm.
- Total weight 0.3 kg.
- Supplied in a carton.

### clar N (Vacuum/LED)

#### clar N – head mirror with LED illumination or vacuum lamp.

- Headband version, with built-in battery compartment and integral charge jack, can be optionally fitted with 2 lithium batteries size CR 123 A, 4 batteries size AAA (alkaline) or 4 rechargeable batteries size AAA (NiMH).
- Headband cleaning made easy by inner, removable and washable padding.
- Particularly comfortable with balanced, infinitely adjustable headband.
- Use without any hindering cables.
- Mirror Ø 55 mm and Ø 100 mm, infinitely adjustable and fixable.
- On/off switch on the battery compartment.
- Charge jack for plug-in charger integrated in battery compartment.
- White LED, 1 Watt with 5.500° Kelvin, 6V or with 6V vacuum lamp.
- 7.300 – 35.000 Lux at 300 mm, depending on illumination version
- 4.000 – 21.000 Lux at 500 mm, depending on illumination version
- Adjustable lamp holder on the mirror to focus the circle of light.
- Battery operation approx. 30 hrs. with lithium batteries, approx. 90 min. with alkaline batteries Type AAA.
- Rechargeable battery operation: approx. 90 min. with fully charged battery.
- Weight with mirror 55 mm: 209 g without batteries, 242 g with CR 123 A lithium batteries, 258 g with rechargeable AAA batteries.
- Weight with mirror 100 mm: 267 g without batteries, 300 g with CR 123 A lithium batteries, 316 g with rechargeable AAA batteries.
- Delivery in carton.

### ri-clar® (halogen)

#### Perfect light reflection with the ri-clar® mirror.

- Models available with mirror diameter 55 mm or diameter 100 mm.
- With 7V halogen lamp for bright, natural light.
- Switch on cable for regulating light intensity in two steps.
- Extremely strong metal ball-and-socket joint permits the mirror to be turned in all directions.
- Manual positioning of the lamp holder for individual adjustment of the reflected light.
- Special fixture located on the side of the forehead for practical and fast exchange of the light heads.
- Simple adjustment and locking of the headband using the rotating aluminum knob.
- Soft, exchangeable foam protection in the forehead area for comfortable wearing of the headband.
- Electrical cable runs along the side of the headband to permit unhindered use.
- Rapid plug-in connector for attaching the headlight cable to the transformer ensures fast implementation.
- Available with 230 V or 120 V transformers.
- Supplied in a sturdy plastic case.

### ri-clar® CL (coldlight)

#### Powerful ENT head light with reflector.

- Precision optics with micro-focus and ball joint for a homogenous and precisely adjustable light point according to size and direction.
- Glass fibre cable for cold-light supply with centred cable duct for undisturbed working.
- Ergonomic headband with head bow and removable leather padding for comfortable use even over several hours, infinite height and width adjustment with central wheels adjustable with one hand.
- Optimum light output by choice of 3 different projectors:
  - HL 150 (150 W halogen)
  - XL 24 (24 W xenon)
  - XL 100 (100 W xenon).
- Fast push connector for connecting glass fibre cable and light source.
- Adapter for projectors of other quality manufacturers.
- Light exit: focusable, double lens system, ball joint.
- Glass fibre cable: length 220 cm, diameter 3.5 mm.
- Total weight 0.3 kg.
- Supplied in a carton.

### ri-magic®

#### Flexible illumination in office practices and hospitals.

- Focused, colour neutral halogen lamp 8V, ideal for examinations and small surgeries.
- The focusing device at the front of the light head allows exact adjustment of the size of the illuminated examination area.
- Heat sink for patient-friendly use.
- Highly flexible, long special arm can be bent in all directions, the light remains fixed in the chosen position.
- Simple to use handle with rheostat for infinitely variable regulation of light intensity.
- 230V or 120V transformers available.
- ri-magic® is available as a desk, wall or floor model
- The floor model is supplied with a particularly steady, heavy and anti-static stand. Two of the five casters can be locked. The wall and desk models can be quickly and easily mounted on all conventional rails using the universal clamp.

### ziegler

#### Head mirror with 90 mm or 100 mm diameters.

- Optional models available with mirror diameter 90 mm or diameter 100 mm.
- Sturdy, easy to assemble ball-and-socket joint to turn the mirror in different directions.
- Simple adjustment to the size of the head using the metal buckle on the side of the headband.
- Supplied in a cardboard box.


Head mirrors, headlights,  
examination lights

## Laryngoscopes

### Riester - Laryngoscopes with unique quality features.


- Exceptional longevity proven by test reports.
- Laryngoscopes available with integrated F.O. illumination, removable F.O. illumination or standard-illumination.
- Excellent XL 3.5 V or 2.5 V xenon illumination and fiber optics or with 2.7V vacuum illumination.
- Macintosh und Miller blades with innovative blade design for optimum visibility of the epiglottis and vocal cords. Simple insertion of the endotracheal tube.
- Large choice of power sources:
  - manoeuvrable and sturdy battery handles and
  - charging bases.
- Wide selection of laryngoscope sets for paediatrics and adults with particularly convincing cost effectiveness.
- All handles and blades are supplied singly and each can be individually combined within the fiber optics or standard model ranges.
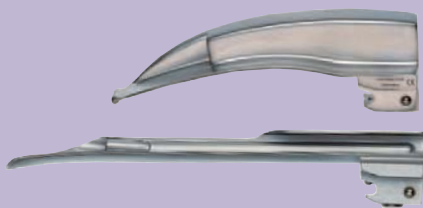
### ri-integral F.O.

### ri-integral flex F.O.
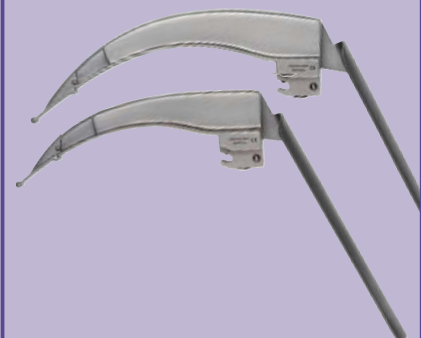
### ri-modul F.O.


### ri-standard

**The results\* prove the excellent light intensity and durability of Riester F.O. blades and their suitability for long-term continuous use.**

\* Test results after 1200 autoclaving cycles (134° C for 5 minutes with the Melag Vacuklav 31). Measurements were performed using a 2.5V XL Xenon lamp at constant light output at a distance of 40 mm from the distal end of the blade (size 3).
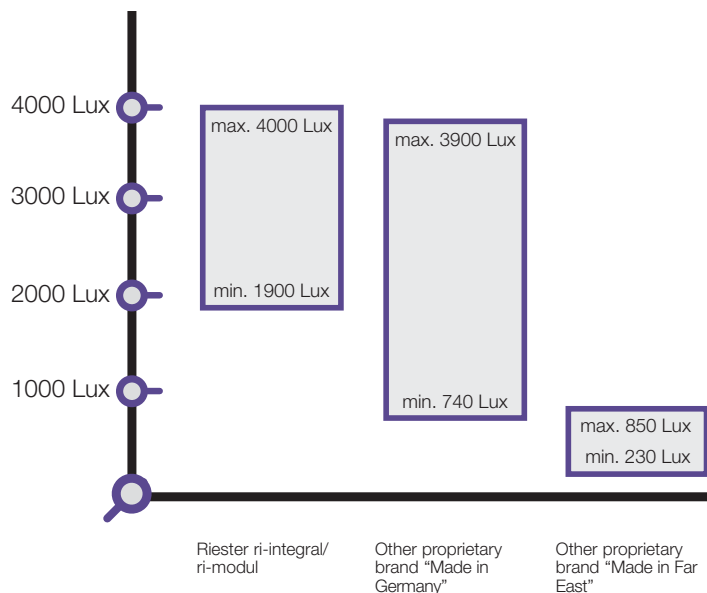
## graeve / cusco

**Riester quality specula. Reliable and durable.**

graeve vaginal speculum

- External polished finish, matt inner surface, available in three sizes (82 mm x 21 mm, 109 mm x 35 mm, 130 mm x 37 mm).

cusco vaginal speculum

- External polished finish, matt inner surface, available in four sizes (for virgins 81 mm x 18 mm, 76 mm x 31.5 mm, 90 mm x 35 mm, 102 mm x 37 mm).
- No reflection thanks to matt inner surface.
- Made of stainless steel type 18/8.
- Can be sterilised up to 134 °C in autoclaves.
- Supplied in a cardboard box.

## speculight

**speculight:- the optimal illumination for all vaginal specula.**

- Battery handle C, chromium-plated.
- For 2 batteries type C.
- 85 cm long cable with fixing clip and 2.5V standard illumination (vacuum).
- Light intensity can be regulated at discretion using the rheostat.
- Supplied in a cardboard box.

## salpingograph

**Tried and tested quality instrument from Riester for the diagnosis of sterility and hysterosalpingography according to Professor Dr. Schultze.**

- For the x-ray contrast imaging of the uterine cavity and the fallopian tubes. Simple test for tubal patency.
- Quality manometer with range 0 up to 200 mm Hg displayed in graduations of 10 mm Hg.
- Highest quality fine mechanical finish with chromium-plated metal parts.
- The instrument can also be read in the dark thanks to the luminescent substance used on the scale.
- Pressure created with a 20 ml syringe via a Luer-lock connector.
- Includes three chromium-plated uterus sealing cones in small, medium and large.
- Holding device for two uterus fixation forceps.
- Complete in high-quality calico case with grey velvet-look inserts.

## Percussion hammers

**Complete range of quality percussion hammers.**

- Reliable and sturdy.
- All hammers supplied in poly-bags.

## Laryngoscopes

## Gynaecological instruments

## Percussion hammers

## Tuning forks

Large range of quality tuning forks and tuning fork sets for ear specialists and neurologists. Selected materials and first-class finish ensure maximum functionality and reliability.

- Highest quality materials: nickel steel or aluminium.
- Large selection of different frequencies: from 64 Hz to 4096 Hz.
- Tuning forks with and without clamps.
- Tuning forks with and without base for different requirements in everyday hospital and office practice.

## Tuning fork sets

The tuning fork sets are well thought-out combinations in an attractive wooden or plastic case to perform full hearing tests. The most important frequencies are all conveniently to hand.

- Set I 8 tuning forks made of steel in wooden case.
- Set II 5 tuning forks made of steel in plastic case.
- Set III tuning forks made of aluminium in plastic case.

## ri-clip®

Fast and simple. The perfect vein compressor for blood extraction.

- Simple blood stasis by tightening the strap.
- Simple release at the clip.
- Avoids unpleasant catching of small hairs.
- Attractive design.
- Colour matching other Riester products.
- Hygienic, can be disinfected with conventional disinfectants or autoclaved at 120 °C/1 bar.
- Available with latex and latex-free!

## komprimeter

Classic device to restrict blood circulation during amputations.

- High precision, durable manometer for accurate pressure monitoring up to 700 mm Hg.
- Sturdy metal pump to produce pressure.
- Air can be released from the cuff via finely regulated, non-ageing air-release valve on the manometer.
- Includes upper arm cuff 57 cm x 9 cm and thigh cuff 96 cm x 13 cm.
- Optional extra: child's cuff 28 cm x 6 cm.
- All cuffs made are made of material which can be washed and disinfected.
- Latex-free.

## spirotest

Exact measurement of lung capacity without the use of water.

- Casing in blue plastic design.
- Measurement range from 1,000 cm³ up to 7,000 cm³.
- Supplied complete in a cardboard box with 50 disposable plastic mouthpieces.
- Suitable for fitness examinations and health checks in the physicians's office or at home.
- Replacement parts: 50 disposable plastic mouthpieces, 1,000 disposable cardboard mouthpieces.

## dynatest®

Simple training and measurement of strength in the hands. For rheumatic disorders or sport injuries.

- Desk device with plastic casing.
- Easy to read scale from 0 to 1.0 bar.
- Black soft plastic bulb.
- Chromium-plated pressure release valve.
- With three yellow coloured soft balls for patients' treatment.
- Supplied complete in a cardboard box.

## metpak®

### For optimum pressure infusions.

- Pressure cuffs with manometer for producing pressure.
- For solutions or blood bags from 500 ml, 1,000 ml, 3,000 ml or 5,000 ml.
- Sturdy cotton cuff in blue, with suspension loop for infusion stands, washable up to 60 °C.
- Precise reading of the level in the bag thanks to a transparent net on the cuff.
- Chromium-plated manometer Ø 49 mm with easy to read aluminium scale up to 300 mm Hg.
- Specially hardened copper-beryllium membrane.
- No zero-point setting required.
- Latex pad resistant to pressure.
- Latex-free inflation bulb.
- metpak® is also available latex-free.

## Veterinary instruments

### Veterinary E.N.T./ophthalmic instrument sets.

- Selection of 2.5V HL halogen and 3.5V XL xenon illumination technology.
- all-metal ri-scope® operation otoscope for a long service life, open design for unhindered use during examination and operation.
- Swivel lens with 2.5-fold magnification.
- ri-scope® ophthalmoscope with glass-fibre reinforced casing, particularly sturdy, light-weight and durable.
- Focussing wheel with corrective lenses from 0 to +40 and 0 to -35 diopters.
- Apertures: Small circle, Large circle, Semi circle, Fixation star, Slit and Grid.
- Filters: Red open filter, Blue filter, Polarisation filter

## med-kit I

### med-kit I

- Sphygmomanometer precisa® N ABS-plastic/metal or alternatively with sphygmomanometer ri-san® (slate-grey, blue, saffron or green).
- Stethoscope duplex® aluminium with replacement parts.
- Diagnostic penlight fortelux® N.
- Percussion hammer buck with replacement parts.

## med-kit II

### med-kit II

- Sphygmomanometer precisa® N ABS-plastic/metal.
- Otoscope uni® I with handle C.
- Stethoscope duplex® aluminium with replacement parts.
- Diagnostic penlight fortelux® N.
- Percussion hammer buck with replacement parts.

## med-kit III

### med-kit III

- Sphygmomanometer precisa® N ABS-plastic/metal.
- Otoscope ri-scope® F.O..
- Ophthalmoscope ri-scope®.
- Battery handle type C.
- Stethoscope duplex® de luxe stainless steel with replacement parts.
- Diagnostic penlight fortelux® N.
- Percussion hammer buck with replacement parts.
- Tuning fork C128 aluminium with clamps.

## Tuning forks

## Products for blood stasis

## Pulmonary pressure measuring devices

## Dynamometers

## Pressure infusion instruments

## Veterinary instruments

## Diagnostic Sets

This brochure shows only a selection of the entire range of Riester medical instruments. We will of course be pleased to send you our fully comprehensive equipment catalog on request. As with all Riester publicity material and sales aids, our comprehensive Riester catalog is also available in German, English, French, Spanish and Russian. All Riester products are supplied with the CE mark. Our general terms and conditions of business apply to the complete catalog.

**Order numbers for the complete catalog:**

German 51231-50 | English 51232-50 | French 51234-50 | Spanish 51236-50 | Russian 51238-50

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